




INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

<p>(51) International Patent Classification<sup>5</sup> : <b>G06F 15/42</b></p>	<p><b>AI</b></p>	<p>(11) International Publication Number: <b>WO 94/11838</b> (43) International Publication Date: 26 May 1994 (26.05.94)</p>
<p>(21) International Application Number: PCT/US93/10525 (22) International Filing Date: 2 November 1993 (02.11.93) (30) Priority data: 07/973,253 6 November 1992 (06.11.92) US (71) Applicant: ABBOTT LABORATORIES [US/US]; CHAD 0377/AP6D-2, One Abbott Park Road, Abbott Park, IL 60064-3500 (US). (72) Inventor: KILLIAN, William, R. ; 939 Wilshire Drive, Libertyville, IL 60048 (US). (74) Agents: GORMAN, Edward, H., Jr. et al.; Abbott Laboratories, CHAD 0377/AP6D-2, One Abbott Park Road, Abbott Park, IL 60064-3500 (US).</p>		<p>(81) Designated States: AU, CA, JP, KR, European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE).  <b>Published</b> <i>With international search report. Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i></p>

(54) Title: PROCESS CONTROL SYSTEM FOR BIOLOGICAL FLUID TESTING

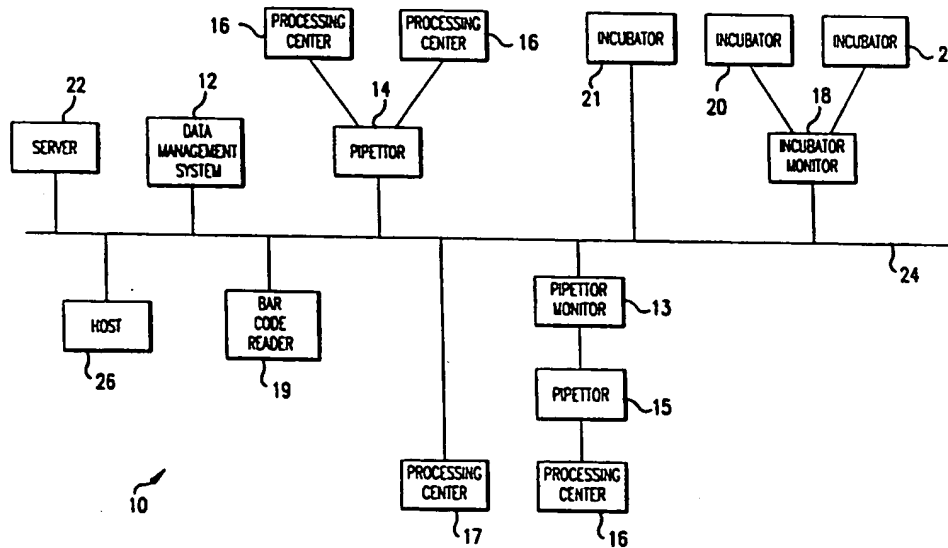

(57) Abstract

A system (10) for automatically testing biological fluids includes a network database (26) which holds test protocols and test results for identified samples. The system (10) tracks the samples by their identifier through the system to testing equipment (16) which is also coupled to the network (24). The system (10) records the tests performed and the times they were performed and monitors any reagents used in the test being run. In addition, the system (10) monitors any incubation of the samples and any results indicated by the test equipment (16). Information on the tests run on a particular sample, the reagents used and the results produced by the tests, are maintained a database and may be formatted into a report.

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Mauritania
AU	Australia	GE	Georgia	MW	Malawi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BC	Bulgaria	IE	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Belarus	KE	Kenya	RO	Romania
CA	Canada	KG	Kyrgyzstan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic of Korea	SD	Sudan
CC	Congo	KR	Republic of Korea	SE	Sweden
CH	Switzerland	KZ	Kazakhstan	SI	Slovenia
CI	Côte d'Ivoire	LI	Liechtenstein	SK	Slovakia
CM	Cameroon	LK	Sri Lanka	SN	Senegal
CN	China	LU	Luxembourg	TD	Chad
CS	Czechoslovakia	LV	Latvia	TC	Togo
CZ	Czech Republic	MC	Monaco	TJ	Tajikistan
DE	Germany	MD	Republic of Moldova	TT	Trinidad and Tobago
DK	Denmark	MG	Madagascar	UA	Ukraine
ES	Spain	ML	Mali	US	United States of America
FI	Finland	MN	Mongolia	UZ	Uzbekistan
FR	France			VN	Viet Nam
GA	Gabon				

- 1 -

## PROCESS CONTROL SYSTEM FOR BIOLOGICAL FLUID TESTING

FIELD OF THE INVENTION

5 The present invention relates to a process control system. More specifically, the present invention relates to a process control system used to document, coordinate, control and verify the testing of biological fluid samples, such as blood, which can be performed by automated test equipment.

BACKGROUND OF THE INVENTION

10 There are many reasons for testing biological fluids. Blood screening of donors at blood banks is a large and still growing market requiring accurate and

- 2 -

reliable testing. Donation of blood typically involves the screening for various diseases and/or foreign bodies which may be present in the blood being donated. In recent years, additional screening tests, such as HTLV-I and HCV have been mandated by law, increasing the number and complexity of tests which each blood sample must undergo. With each test, additional FDA minimum standards have also been set. Other biological fluids such as urine and saliva tests are routinely tested and/or screened for substances, diseases and chemicals.

Over the years, as both the amount of blood being donated and the number of required tests have increased, a degree of automation has occurred. Systems such as the Commander System<sup>TM</sup> by Abbott Laboratories, have partially automated many of the steps involved in blood testing.

The Commander System can be divided into four main parts: the Flexible Pipetting Center (FPC), the Parallel Processing Center (PPC), the Dynamic Incubator (DI), and the Data Management System (DMS). Under this system, a donor's blood sample is initially bar coded, linking the blood sample to a database of information about the patient, including the test(s) to be run on the blood specimen. This information is stored in a local database with the same information residing in the DMS database. The bar coded samples are then taken over to the FPC for pipetting into test wells in which the testing of the blood specimens will take place. The bar coded blood samples, stored in test tubes, have their exact position in a

- 3 -

holding tray entered into a database. This position information is used by the FPC in pipetting the samples.

5 The test wells are typically located in test trays. Each tray may contain 60 or more test wells into which blood specimens can be pipetted for testing. The FPC, at the same time that the pipetting of blood specimens takes place, can add other materials such as a diluent, as may be required by a particular test. Each tray is also bar coded. This bar code is read as soon as the tray is placed into its appropriate receptacle in the FPC. The bar code information on the tray, combined with the bar code information on each blood sample, allows accurate data gathering to track the blood sample through the testing process. In the example of the Commander System, up to 15 four trays (and thus four tests at one test per tray) can be placed in the FPC at one time. The FPC has a maximum throughput of 900 tests per hour, highlighting the importance of accurate coordination and tracking.

20 Once a blood sample has been pipetted into the test wells and any required diluent has been added, the tray containing the test well moves either to the Dynamic Incubator or to the Parallel Processing Center, depending upon the tests being run on the particular blood specimens. In many cases, beads coated with an antigen or antibody, 25 required for the test, are also added to the test wells to aid in the detection of various substances.

At the Parallel Processing Center, the bar code of the tray is again read, to insure tracking of each

- 4 -

individual blood sample contained in the test wells. This provides reliability in linking test results to each individual blood sample. A connection, such as a serial data connection, links the FPC to the PPC to track the blood sample through the testing which takes place at the Parallel Processing Center. In the Parallel Processing Center, washing of samples, as well as addition of reagents takes place according to a preprogrammed protocol or "recipe" for each specific test which the samples in each tray undergo. Should there be any problem with reading the bar code, an operator may physically enter the tray identification number to ensure accurate tracking. The PPC indicates the completion of the processing pass by way of an audible tone. The PPC also reads the test results by way of a spectrophotometer. The test results, which are linked to a particular blood sample, by way of the bar code, are then ported over to the DMS via a connection, such as a serial data connection.

The Dynamic Incubator provides incubation and/or agitation for samples undergoing tests which require these steps. The DI of the Commander System can incubate up to eight trays at various possible preset temperatures (although only one temperature at a time may be set). The DI also provides for entering specific times for each individual tray placed in the incubator. When the entered time has elapsed for each individual tray, a signal, such as an audible alarm combined with a display, will indicate which tray has completed its incubation/agitation.

- 5 -

The Data Management System takes the information provided to it from the FPC and PPC to provide a report of the test results for each blood sample.

5 Although blood testing systems, such as the  
Commander System, have provided a degree of reliability,  
accuracy and consistency in the blood testing process, the  
increase in amount of blood donated, the number of tests  
which can be run for each donor, and the ever increasing  
10 government standards, including the Good Manufacturing  
Practices (GMP), make even better safeguards and better  
control and documentation of the entire process more and  
more desirable. Further, the increase in the number of  
tests which can be performed or are required to be  
15 performed, means additional chemicals and reagents which  
must be closely monitored for factors such as expiration  
dates. This raises the possibility of operator error in  
selecting, tracking and operating the tests. Lastly, the  
ability to monitor the integrity of each test and the  
ability to determine whether multiple tests (such as  
20 retests to provide statistical accuracy) are required are  
also desired features.

#### SUMMARY OF THE INVENTION

25 The present invention is embodied in a testing  
system for biological fluids which includes means for  
transferring the fluids from a test tube or holding  
container to a test site, such as a test cell in a testing

- 6 -

tray. Means are also provided for testing the biological fluid samples for a number of desired tests. Control means, such as a computer, coordinate the transfer and testing of the biological fluids, as well as tracking these operations. Lastly, connection means are provided, allowing communication among all of the devices.

In an exemplary embodiment of the present invention, a Pipetting System is used to deliver blood samples, for example, to a test well in a test tray which is the test site. Either or both a Processing Center and an Incubator (under control of an Incubator Monitor) are used as required by the specific test which the blood sample is undergoing. A computer, such as a personal computer, coordinates the testing, tracking bar code information, directing equipment and human operators to ensure that the proper test procedures are carried out and the proper chemicals and reagents are used. At the same time, the computer documents whether each test has been carried out as specified, including whether the test equipment has been calibrated and whether the chemicals and reagents used are not beyond their expiration dates. A network, such as an Ethernet<sup>TM</sup> network, connects the various pieces of equipment to the computer for effective communication among each piece of equipment and the control computer.

BRIEF DESCRIPTION OF THE FIGURES



- 7 -

The present invention is described by way of nonlimiting example with reference to the following figures in which:

5 Figure 1 is a system diagram of an exemplary testing system of the present invention;

Figure 2 is a program flow diagram which illustrates the operation of the TIP/TRAY LOT NUMBER SPECIFICATION routine;

10 Figure 3 is a program flow diagram which illustrates the operation of the MASTER LOT INFORMATION routine;

Figure 4 is a program flow diagram which illustrates the operation of the VIEW EXPIRED MASTER LOT INFORMATION routine;

15 Figure 5 is a program flow diagram which illustrates the operation of the MASTER LOT INFORMATION routine;

20 Figure 6 is a program flow diagram which illustrates the operation of the assay section of the PIPETTING RUN SETUP routine;

Figure 7 is a program flow diagram which illustrates the operation of the control section of the PIPETTING RUN SETUP routine;

- 8 -

Figure 8 is a program flow diagram which illustrates the operation of the destination object specification section of the PIPETTING RUN SETUP routine;

5 Figure 9 is a program flow diagram which illustrates the operation of the diluent and control placement section of the PIPETTING RUN SETUP routine;

10 Figure 10 is a program flow diagram which illustrates the operation of the volume verification of diluent and control section of the PIPETTING RUN SETUP routine;

Figure 11 is a program flow diagram which illustrates the operation of the PIPETTING routine;

15 Figure 12 is a program flow diagram which illustrates the operation of the MONITOR AND FUNCTION ENTRY routine;

Figures 13a and 13b are program flow diagrams which illustrate the operation of the BEGIN TRAY INCUBATION routine;

20 Figures 14a and 14b are program flow diagrams which illustrate the operation of the TRAY INCUBATION COMPLETE routine;

Figures 15a and 15b are program flow diagrams which illustrate the operation of the ALARM STATE routine;

- 9 -

Figure 16 is a program flow diagram which illustrates the operation of the PRINT INCUBATOR DATA routine;

5 Figure 17 is a program flow diagram which illustrates the operation of the PRINT INCUBATOR MONITOR ERROR LOG routine;

Figure 18 is a program flow diagram which illustrates the operation of the TIME OUT OF COOLER routine;

10 Figure 19 is a program flow diagram which illustrates the operation of the COVER SEAL LOT NUMBER ENTRY routine;

15 Figures 20a and 20b are program flow diagrams which illustrate the operation of the BEAD ADDITION routine;

Figures 21a and 21b are program flow diagrams which illustrate the operation of the CONJUGATE RECONSTITUTION routine;

20 Figures 22a and 22b are program flow diagrams which illustrate the operation of the OPD PREPARATION routine;

Figure 23 is a flow chart diagram which illustrates a three-pass test which may be performed using the system shown in Figure 1;

- 10 -

Figure 24 is a program flow diagram which of the PREPARE REAGENTS routine shown in Figure 23;

5 Figure 25 is a program flow diagram of the CHECK VALID TRAY routine which is part of the PROCESS TRAY routines shown in Figure 23;

Figure 26 is a program flow diagram of the APPROVE TECHNICIAN I.D. routine which is part of the PROCESS TRAY routines shown in Figure 23;

10 Figures 27a and 27b are program flow diagrams which illustrate the remaining steps in each of the PROCESS TRAY routines shown in Figure 23;

Figure 28 is a flow diagram which illustrates a collation of results by the DMS;

15 Figure 29 is a flow diagram which illustrates exemplary report options;

Figure 30 is a flow diagram which illustrates exemplary batch record and statistical reports;

20 Figure 31 is a flow diagram which illustrates a test which may be performed using the system shown in Figure 1.

DETAILED DESCRIPTION OF THE INVENTION

- 11 -

There is shown in Figure 1 an exemplary system for testing biological fluids or Process Control Network 10. Testing system 10 comprises a Data Management System (DMS) 12, a Pipettor 14, a Processing Center (PC) 16, an Incubator Monitor (IM) 18 and an Incubator (I) 20. There are multiple units of Processing Center 16 and Incubator 20 shown to illustrate that more than one of these units (as well as the others) can be connected in system 10. There is also shown in figure 1 a Pipettor 15 which is connected to a Pipettor Monitor 13 to illustrate that the "intelligence" for a pipettor can be separate from the pipettor. Similarly, Incubator 21 is shown as being directly connected to network 24 to illustrate that the intelligence or front end processing can be a part of the incubator. Other configurations, not shown, are possible depending upon the needs for system 10 by the user.

Bar Code Reader (BCR) 19 is shown for reading bar code information and sending the information out onto network 24. A bar code reader can also be connected to a particular device for entering bar code information at particular locations on the system and/or times in the process.

Also shown is a host computer 26, server 22 and a network 24. Host computer 26 may be a bulk storage computer or mainframe for holding a large database of patient information and results from past tests. Network 24 may consist of an Ethernet network for example, which connects each of the devices to server 22. A network such as an Ethernet requires network cabling, network cards, a

- 12 -

network operating system and assorted hardware as is understood by those skilled in the art. Numerous other types of networks can be used to allow communication between the various devices, making up the present invention. Similarly, server 22 may consist of any computer suitable for operating the network being used. A personal computer using a 386, 486 or similar microprocessor chip can be used as server 24. RISC workstations or mini or mainframe computers can also be used to operate the software and control the network which connects all of the devices of the present invention to provide communication.

The components of the present invention are described by way of figure 1 using those components of the Commander System (incorporated by reference, herein) which are applicable to the present invention. This is in no way meant as a limitation of the equipment which can be used. For example, pipetting systems other than the Flexible Pipetting Center can be used, controlled and tracked in accordance with the present invention. Similarly, incubators other than the Dynamic Incubator may be used. The following are Abbott Laboratories List Numbers for components which can be used in an exemplary embodiment of the present invention:

25	<u>Item</u>	<u>List #</u>
	Flexible Pipetting Center 14	3A46-01
	Personal Computer 12	3A46-70
	Control Software (not shown)	3A46-87
	Bar Code Reader 19	3A46-28
30	Dynamic Incubator 20	6210-01
	Parallel Processing Center 16	6208-01

DMS 12 is used to enter and coordinate information regarding the tests to be run for a given  
5 sample. The DMS relays this information via network 24 to the devices which carry out the testing steps. DMS 12 operates a database, such as Oracle<sup>TM</sup> Version 6.X to track, store and answer system requests for information.

DMS 12 evaluates the information for the samples  
10 (or requests to test samples) and translates this into the actual tests to be carried out. If a sample is to be tested for diagnostic purposes, the appropriate tests would be ordered, which may differ from the tests ordered for a blood donor screening. For instance, a blood sample from  
15 an HIV patient would not be subjected to routine screening tests, but rather to tests related to the HIV condition which have been predetermined and entered into system 10. Conversely, a normal blood donor sample would be given a preselected group of tests regardless of the condition of  
20 the donor, because it is a screening process.

Reagents are used in testing a biological fluid sample. Reagents are logged into a reagent library containing information on each reagent lot, including  
expiration dates, lot numbers and amounts of reagent  
25 received. The reagent library is a database of information stored on a computer. The database can be set up using commercial products, such as Oracle<sup>TM</sup>. It is also possible to provide any particular information necessary for use of

- 14 -

the reagents, such as optimum storage and operation  
temperatures, and parameters to use to normalize the  
results of tests performed using the reagents to results  
obtained from an ideal reagent. The amount of information  
5 logged for any reagent lot would be selected by the user  
and may depend upon factors such as the tests being run  
and/or government regulations.

In an exemplary embodiment of the present  
invention, the reagent information accompanies each lot of  
10 reagents. Each bottle of reagent has a bar code label  
affixed to the bottle, linking the bottle to the master lot  
information contained on the sheet which accompanies each  
kit of reagents. The lot of reagents and the master lot  
information sheet comprises a reagent kit. As a reagent  
15 kit is received, the bar code information from the master  
lot information sheet is scanned into the network via a bar  
code reader, such as BCR 19. This information is then  
stored on the network, in a local memory storage device,  
such as a hard disk drive, where it can be made available  
20 to all equipment via network 24.

For example, the reagent kit information can be  
read into the system using a bar code reader attached to  
a Pipettor 14 and stored on the hard drive in the control  
computer of Pipettor 14. Another embodiment may have the  
25 information stored at a central location such as server  
22. These and other variations will be understood by those  
skilled in the art.



- 15 -

A bar code reader may be separate from or attached to a Pipettor 14. As the system of the present invention is connected via a network, the bar code information can be read in from anywhere on network 24 and then stored at its proper location.

Once the master lot information is stored in system 10, reagents can be used until an expiration date, for example, is noted by system 10. At this point, the remaining reagents in the lot, if any, would be removed and a new reagent kit brought in. If system 10 detects a reagent bar code which is not associated with an approved reagent, the operator will be alerted and the processing halted. At this point, another approved library lot must be selected for use or a new acceptable lot must be loaded into the database such that it can be subsequently used for processing.

Automatic test ordering can take place with the present invention through host computer 26. Host computer 26 stores the original donor information. It is not necessary and often not advisable for a testing facility using the present invention to possess all donor information which may include, for example, the name of the patient. The present invention only requires information necessary to carry out the testing of biological fluid samples such as blood. Host computer 26 can access DMS 12 directly through network 24, requesting specific tests for a particular blood sample. Host computer 26 can also be connected to DMS 12 via a serial connection (not shown), as opposed to the network connection. As DMS 12 initiates

- 16 -

and coordinates the testing process, there would be no need for the individual pieces of test equipment to access host computer 26 directly or vice versa. It is nevertheless possible for all equipment to access host 26 directly. In  
5 any event, the particular tests designated for a given blood sample would be communicated from host computer 26 to DMS 12 for storage and eventual execution by system 10.

Pipettor 14 is used to pipet an aliquot of sample from a test tube into a test site. Aliquots of  
10 sample are typically in the range of less than 300 micro liters, but amounts may vary, depending upon the test involved. In this example, the test sites are the test wells of a test tray which fit into receiving portions of Pipettor 14, Processing Center 16 and Incubator 20. In  
15 this example, each test tray has approximately 60 test wells and thus 60 samples can undergo the test being applied to the tray.

For this exemplary embodiment, Pipettor 14 as was used in the Commander System is being described. This  
20 pipetting system uses a single pipet with a removable pipet tip. It is well within the scope of this disclosure that other pipetting systems can be used within this invention. Other pipetting systems may use multiple pipets which can pipet up to 12 samples at one time. Many bulk applications  
25 use a multiple pipet system for delivery of sample to multiple test sites.

Pipettor 14 polls DMS 12 to obtain information necessary to pipet an aliquot of sample into the proper

- 17 -

number of test sites. For example, if four tests are desired, then four trays would be placed in Pipettor 14 (where each tray of test wells undergoes one test protocol). DMS 12 communicates with Pipettor 14 over network 24, providing the test information for each sample to Pipettor 14. This can be done for each individual sample or for blocks of samples. Pipettor 14 pipets an aliquot or portion of the sample into one test well per tray to provide four tests. Should fewer tests be necessary for a particular sample, Pipettor 14 would only pipet according to the number of tests necessary. Other samples may require fewer or more tests. By having this information in DMS 12 and, having Pipettor 14 in communication with DMS 12, appropriate tests are ordered for each sample. An aliquot of a given sample is then pipetted to the proper test site in the proper tray automatically. There is thus, no requirement of operator intervention to allocate the sample to the proper test trays nearly eliminating errors caused by operator intervention. This is the implementation of the sample profile originally downloaded from the host or entered into the system.

Pipetting can begin at any time in the process subsequent to downloading sample information from host 26 to DMS 12. The reagents are bar coded and entered into the system as described above. Although the present embodiment is described using bar codes, other methods of encoding and labelling could be used to provide easy entry and tracking of samples, trays and reagents. Each Pipettor 14 checks to make sure that proper, unexpired reagents (such as

- 18 -

calibrators, controls, or other reagents) are present before it begins pipetting a sample into the designated test sites of the trays. If there are invalid reagents due to situations such as an expired reagent lot, no pipetting will occur. Each reagent volume will be checked prior to pipetting to ensure adequate volume is present for the batch.

Further, in running the test, it is a common practice to use control specimens along with the sample specimens. The control specimens are used to ensure that the reagents and/or other testing facilities are working properly to identify the substances for which the samples are being tested. It is typical to run a control at the beginning of a group of trays being subjected to the same test for a particular substance or disease. Thus, Pipettor 14 also makes sure that there is sufficient control specimen, before a test is commenced.

The test tubes of samples are loaded onto Pipettor 14 and an automatic bar code reader (not shown) scans the sample tubes. It then polls DMS 12 as to which tests are appropriate for the given samples. This information, combined with the information concerning the reagents which are loaded on Pipettor 14 at the present time, determine the tests which the samples will be prepared for on Pipettor 14 at that time. For testing situations which require multiple runs to ensure proper statistical results, system 10 would be advised of this through network 24 with information from DMS 12 and

- 19 -

instruct Pipettor 14 to pipet the aliquots of samples into multiple test sites as appropriate.

5 In essence, the present invention removes the need for personnel to manually check reagent expiration dates, presence of proper reagents in the system, proper reagent volume and other factors. Further, the ability to run multiple replicates of a given sample within a batch alleviates the need to sort sample tubes into those requiring multiple tests and those that do not. The use of  
10 the bar codes combined with the communication between master database and testing equipment allows preparation of a variety of aliquots of samples requiring different tests.

Processing reagents for PC 16 may take the form of conjugate, acid, OPD (o-Phenylendiamine, a color tag),  
15 or coated beads, for example. As each reagent is loaded onto a Processing Center 16, the bar code information for the reagent is scanned using a bar code scanner attached to each Processing Center 16 or Incubator Monitor 18.

The reagent library is stored on system 10. It  
20 can be stored at a central location, such as server 22 or DMS 12 or it can be distributed at multiple locations such as Pipettors 14 and Incubator Monitors 18. As the information is entered, it is also checked to make sure that the reagents have not expired. This checking will  
25 take place prior to use on each batch, although a user could designate different check times depending on the tests being run.

- 20 -

As is shown in Figure 1, there are multiple Processing Center 16 units which can be connected to each Pipettor 14 and/or network 24. In the exemplary embodiment shown in Figure 1, only Pipettor 14 is connected to the  
5 network. It is wholly within the scope of this invention and disclosure that any and all devices can be directly coupled to the network for communication between the particular devices and DMS 12.

In this exemplary embodiment, Processing Center  
10 units 16 are connected to network 24 via Pipettor 14, reducing the need for network communication equipment in the Processing Center units 16. A Processing Center 17 is, nevertheless, illustrated to show that such a configuration can be used. It is also possible for any Processing Center  
15 16 to poll network 24 for necessary information stored in any Pipettor 14. This allows constant and accurate tracking, as well as reliable testing at any point in the testing process without the need for operator input of testing information. Each device can request and receive  
20 the information it needs to carry out its functions.

With pipetting complete, the tray of test sites is transported to a Processing Center 16 for processing. Trays may be transported by operator or mechanical means, such as a conveyor belt. This applies to the transport of  
25 trays (and thus test sites) between any devices in system 10. As the tray enters Processing Center 16, the bar code on the tray is automatically read and a request is sent out over network 24 for assay information. The pipetting information can also be stored locally at a particular

- 21 -

Pipettor 14 for communication to Processing Center 16 units. Processing Center units 16 use this information to carry out the assays using the reagents loaded into each Processing Center 16.

5           In many instances, it is necessary to provide several passes through a Processing Center 16 to complete the assay protocol steps for a particular test. It is also typical for tests to undergo incubation after each pass through a Processing Center 16.

10           For instance, after a first pass through a Processing Center 16, it may be necessary to incubate the tray of test sites. At this point, the trays would be transported from Processing Center 16 to an Incubator 20 where the tray bar code is scanned (either at Incubator 20  
15 or IM 18) prior to placing the tray in one of the selected tray slots. Incubator 20 uses the network to check the tray status and confirm incubation requirements. In the present embodiment, information is communicated to an Incubator 20 through IM 18. In effect, IM 18 possesses the  
20 "intelligence" and processing capability, similar to the computer controlling each Pipettor 14. The incubation requirements are based on the particular samples being tested and the particular test being run for those samples. The necessary incubation time and temperature would be part  
25 of the profile for a particular test. A IM 18 requests the information linked to the bar code through network 24. The request would obtain the information from a Pipettor 14 or other location where the information is stored.

- 22 -

Once a tray is inserted into an Incubator 20, its position is confirmed by scanning a bar code adjacent to the incubation slot and beginning incubation. An Incubator 20 will automatically be set to the proper incubation time, temperature and agitation of the tray (if necessary) by IM 18. IM 18 will constantly check temperature to confirm that all specifications for the tests have been met.

One of the features of the present invention is the ability to catalog and manage all of the data which is present during the testing of aliquots of samples. Information such as incubator serial number, time of incubation, reagent lot, etc. are now all stored and can be retrieved at any time. For instance, if at a later time it is determined that a particular lot of reagents was actually defective, all testing done with that lot can be immediately identified and flagged as invalid. Also, should testing requirements change, a system wide change can be implemented at any time. This would also allow the flagging of samples done under prior methods which may prove invalid under new specification. Data is available to the system during the testing process. The length of time the information is stored on system 10 depends upon the particular user and capacity of equipment such as hard drives, large volumes of information may be deleted at selected time intervals or archived in an archive storage device (not shown). There is little need for operator intervention as the desired information is already programmed into the system for retrieval and monitoring. Operator input is basically limited to scanning bar codes



- 23 -

where this is not automatically done and entering operator identification codes as steps are performed.

5 Because test results are a part of the database information, testing which requires multiple runs for statistical purposes can be automatically implemented once the first test results are obtained. Once a first positive result is obtained, the existence of this result in the database will cause system 10 to automatically run the necessary additional test for the sample, the next time the  
10 sample is loaded into Pipettor 14.

Once all testing is complete, the information from the testing cycle can be put into a report sorted and restricted as is understood by those skilled in the art. This information can be passed along the network to other  
15 facilities, maintained in a data storage facility, printed out by a printing device, or displayed on a display device.

The following routines are illustrative of the exemplary embodiment and should not be read as solely limiting the operation of the various components of the  
20 present invention.

There is shown in figure 2 a routine which addresses the specification of the tip or tray lot numbers. In step 32, either the tip or tray lot number is selected from the system menu. The exemplary embodiment of the  
25 present invention is menu driven for ease of operation by the operating personnel. Accordingly, many of the functions will be carried out through selection of menu

- 24 -

choices. Implementation of a menu drive selection system as well as other systems, will be understood by those skilled in the art.

5 In step 34, the operator selects the pipettor to be used if more than one pipettor is configured in the system. In step 36, the operator would enter his or her operator identification information. This is for tracking purposes, so that the personnel carrying out the various tasks requiring operator interface can be tracked. In step 10 38, an electronic work sheet (configuration menu) is produced for the operator listing the locations of all disposables and their current lot numbers. Any lot numbers that have changed should be updated on the system. It should be noted that the entry of lot numbers must take 15 place prior to their use in system 10. In step 40, the selections are saved into system 10, updating the database of information which is stored on system 10 at an appropriate storage location.

20 There is shown in figure 3 a routine for viewing the master lot information. In step 44, the operator selects the master lot information choice from the main menu. In step 46, the operator selects the view master lot choice from the master lot information menu. This is a menu nested within the main menu. In step 48, the operator 25 goes to the next nested window wherein the master lot which is to be viewed is selected from a list of all master lots available in the system. In step 50, the operator views the list of all information associated with the desired master lot. Then in step 52, the user may delete the

- 25 -

5 master lot and all information associated with it. This will remove the master lot from the logical area where it could be used by system 10. Deleting of any information will typically be password controlled, at the option of the system supervisors. Master lot information may be deleted if at any time, for example, the quality of the reagents in the master lot becomes questionable.

10 There is shown in figure 4 a routine for viewing expired master lot information. In step 56 the operator selects the master lot information choice from the main menu. In step 58, the operator selects to view the expired master lot choice from the master lot information menu. In step 60, the operator receives a report of all expired master lots in the system. This can be either on screen  
15 (not shown) or by a printout of a printer (not shown) attached on the network. Finally, in step 62 the user is given the option of deleting the expired master lots. Again, for tasks such as this, password control will probably be desired.

20 There is shown in figure 5 a routine for recording master lot information. In step 66, the operator selects master lot information from the main menu. In step 68, the operator selects the record information option from the master lot information menu. In step 70, the operator  
25 enters the password if password protection is enabled. In step 72, the operator enters his or her operator identification code to maintain tracking of users. In step 73, the master lot number from the kit card is entered. This card is enclosed with the reagent kit by the

- 26 -

manufacturer. As previously described, this card will contain pertinent information to the reagent lot being entered. In step 74, a determination is made as to whether the master lot information has already been entered. If so, the current information is displayed and processing moves to step 84. Otherwise, processing continues to step 76 where the test identifier for the master lot is entered from the kit card. The test identifier provides the information for the tests to be run with the reagents from the kit. In step 78, the number of components in the kit is then determined. In step 80 the lot number is entered for each component of the kit. In step 82 the expiration date is entered for each component of the kit. In step 84 a loop back to step 78 is performed until all of the component information from the kit has been entered. In step 86 the operator requests confirmation of the components in the kit to check against the kit listing. In step 88, the master lot information is stored in a database which is accessible via network 24. In step 100, the master lot summary report is printed.

The routine for an exemplary assay selection is shown in Figure 6. The pipetting run setup begins with the assay selection. In step 104, the pipetting setup is selected from the menu. In step 106, the operator's I.D. is entered for the pipetting run, about to begin. In step 108, the sample tube inner diameter is entered. This provides the necessary information for Pipettor 14 which in the exemplary embodiment is the Flexible Pipetting System 14 which uses stepper motors to take the sample fluid and deliver it to the test sites. In step 110, the bar code

- 27 -

mode is entered for reading bar code information. In step 112, a request is made to select the tray map print option, should the operator desire to print out a map of the tray. In step 114, the pipetting type is selected; for example  
5 batch or random mode. In step 116, the selection of the assays to be pipetted during the run is selected. In step 118, a decision is made whether the assays can be pipetted together. If so, processing moves to step 122. In step 120, a decision is made whether assays can be pipetted  
10 together and if not, processing returns to step 116.

In step 122, the number of destinations for each assay is determined. This corresponds to the number of tests which will be run on each fluid sample. In step 124, it is determined whether Pipettor 14 is being run in the  
15 total process control (TPC) mode and whether Pipettor 14 is in the record or verify mode. If so, the master lot number for the assay is entered.

To allow for flexibility, network 24 may be configured for total process control (TPC) OFF/ON. If ON,  
20 the user may then further turn ON/OFF various features as he or she desires. For example, one user may want the system to track and verify reagent status but may not wish to enter and track technician identification.

Furthermore, the user may want to run certain  
25 tests under total process control, as this was defined by configuration described above, but they may want to run other tests without total process control. In order to

- 28 -

provide for this flexibility, the user may define modes for each assay as follows:

Normal Mode - no features of total process control are enabled;

5                   Record Mode - for those features turned ON, system 10 will record and track data as it is entered. However, this data will not be verified. For example, if a reagent lot number is entered, it will be recorded but not verified in the reagent library. The user may also choose  
10 to enter no lot number in which case the system will move on to the next step without alarming the user.

Verify Mode - for those features turned ON, the system will record, track, and verify data entered. If  
15 data is found to be unacceptable, or if no data is entered, the user will be notified with an alarm and be prompted to correct the problem. The user may choose to override the alarm state, in which case the data is then only tracked and is flagged as a deviation.

In step 126, if the master lot number does not  
20 turn up, the operator will be questioned as to whether this is a reagent qualifying run. If so, the information concerning this assay will be flagged to indicate that it is in qualifying mode. This will provide calibration and testing information for the system. In step 128, if the  
25 TPC mode is ON, the verify mode is ON, the qualify mode is OFF, and the master lot number does not exist in the database, the system will provide a warning message to the

- 29 -

user. An option of overriding the warning is provided allowing this step of the assay to be processed in record mode. If the operator chooses not to continue, processing will return to step 116. In step 130, a determination is made if the requested number of destination objects for each assay causes the destination objects not to fit on the platform. In such a case, processing also returns to step 116. In step 132, the platform configuration is saved.

In Figure 7, there is shown the control section routine. In step 134, the controls to be pipetted for each destination object during a run are selected. The controls are used to verify the accuracy of the testing equipment and reagents. In step 136, the control selection is saved.

In Figure 8, the destination object specification is provided. In step 138, a determination is made as to whether the assay has a reagent coated destination object, the system is in TPC ON mode and the assay mode is set to record or verify. Under such conditions, the operator is asked to enter the lot number of the destination object. In step 140, the user enters the actual destination object. In step 142, if the destination object I.D. exists in the system, another message is provided and processing returns to step 140. In step 144, a new destination object I.D. is created in the database. In step 146, the user is asked to place the destination object on the platform and confirm that proper placement has taken place. In step 148, a determination is made as to whether there are more destination objects that need to be placed. If so, processing returns to step 138.

- 30 -

Otherwise, processing continues to the routine shown at Figure 9.

Figure 9 shows the routine for the diluent and control placement. In step 150, a screen is displayed  
5 indicating the location for all reagents and controls on the platform of Pipettor 14. In step 152, a determination is made as to whether the system is in TPC ON mode and the assay mode is set to record or verify. If so, a request of the lot number for the diluent or control bottle is made.  
10 The user may override these modes. In step 154, the I.D. of the control (should a control be placed at this time) is made against the master lot list. In step 156, if the I.D. does not match an entry in the list, the operator is provided with an error message and processing returns to  
15 step 152. In step 158, if the expiration date of the control has been reached, the operator is provided with an error message and processing proceeds to step 168. In step 160, a determination is made as to whether the master lot of the control is consistent with the other controls in use for this assay. If inconsistent, the operator is provided  
20 with an error message and processing jumps to step 168. In step 162, the operator is asked to place a control bottle in the proper selected position. In step 164, the operator is instructed to bar code the position of the control  
25 bottle. If the position differs from the expected position, the operator is provided with an error message and this step is repeated. In step 166, processing jumps back to step 152 and is repeated for all remaining diluents and controls for this pipettor run. If all diluents and  
30 controls have already been assigned positions, processing



- 31 -

proceeds to step 174. In step 168, the operator is given the opportunity to override the error message. If this option is taken, processing returns to the earlier. If the error is not overridden, processing proceeds to step 170.

5 In step 170, the operator is requested to remove the control bottle from the platform. In step 172, processing goes back to step 152. In step 174, the platform configuration is saved and in step 176, a determination is made as to whether the dilutor should be primed.

10 There is shown in Figure 10 a routine for volume verification of diluents and controls. In step 178, a determination is made as to whether the TPC mode is in the off position. If so, processing proceeds to step 214. If not, a new pipet tip is picked up from the rack in step  
15 182. In step 184, the fluid level in the diluent bottle is sensed to determine the amount of diluent which is present. In step 186, the number of wells that could require diluent for this particular assay is computed. In step 188, a warning is provided should the amount of diluent be  
20 insufficient to fill an entire tray as determined in step 186. In step 190, the operator is given the option of replacing the bottle of diluent.

Should the operator choose not to, processing proceeds to step 198. In step 192, the operator is asked  
25 to remove the diluent bottle from the platform. In step 194, if the assay mode is set to verify or record, the operator is asked to enter the lot number for the diluent. In step 196, if the assay mode is set to verify, the lot number on the master lot list is verified. If incorrect,

- 32 -

the operator is given an error message and processing returns to step 192. If correct, the user is instructed to place the bottle of diluent on the platform. In step 198, if more diluents are present, processing returns to step 5 180. Otherwise, processing continues to step 200.

In step 200, a determination is made as to whether the assay used in the control bottle is set to normal mode. If so, processing proceeds to step 212. Otherwise, processing continues to step 202 where a new tip 10 is picked up from the rack. In step 204, the fluid level in the control bottle is sensed to determine the amount of fluid present. In step 206, the tip is dropped into the waste bin for disposal. In step 208, if sufficient fluid is available in the control bottle to pipet all control 15 positions for this control in this assay, processing continues to step 212. In step 210, a determination is made as to whether the fluid available is adequate to pipet all control positions. Pipettor 14 will keep track of the bottle position and the amount of fluid required. Volume 20 verification is continued for the remaining bottles.

In step 212, a determination is made as to whether there is sufficient fluid and whether there are more control bottles to fill. If so, processing goes back to step 200 for the remaining control bottles. In step 25 214, if all bottles have sufficient volume, pipetting begins with the routine set forth in Figure 11.

In step 216 of figure 11, a determination is made as to whether any bottles do not have sufficient

- 33 -

5 volume. If so, the information is displayed to the operator and the pipettor is stopped. In step 218, the operator is given the option of overriding the warning and beginning the pipetting. In step 220, the operator is asked to remove all control bottles for all affected assays from the platform. In step 222, processing returns to step 150 in Figure 9. This ends the pipetting run setup routines.

10 There is shown in Figure 11 a routine for pipetting in random mode. In step 224, Automatic Bar Code Reader (ABC), the additional device attached to the side of Pipettor 14 to hold and read tubes, is indexed to read a sample. If the sample is the last sample of a batch, processing is stopped. If not, processing continues to  
15 step 226 where the bar code for the sample is read.

In step 228, a look up of the bar code takes place for the sample table. In step 230, a determination is made as to whether the sample is to be pipetted to assays currently being processed and a determination is  
20 also made as to the number of replicates for each destination object. If the sample is not to be pipetted, processing returns to step 224. Otherwise, processing continues to step 232 where the correct number of replicates is pipetted to each destination object. The  
25 number of replicates is to be stored with the sample I.D. in the sample database, rather than with the assay protocol or test procedure. The assay protocol is the processing algorithm or "recipe" used for a particular test (assay). In other words, the specific series of steps and data

- 34 -

reduction parameters associated with a particular assay. Processing then returns to step 224 until all pipetting is complete. For certain steps (as those described with respect to figures 2-11) the operator is prompted to enter  
5 his or her identification code. System 10 then determines if the operator is qualified for the task at hand.

Although not shown, a network clock signal can be broadcast to each Pipettor 14 from DMS 12 (step 234, figure 11). Pipettor 14 in turn, updates its internal  
10 clock as well as forward the network clock message to all PC 16 units which are connected to it (step 236, figure 11). In the case of PC 16 units directly connected to network 24, they would receive the network clock message independently.

15 The main subroutine for the operational phase of the system is the monitor and FUNCTION ENTRY routine. This routine, which runs on incubator monitor 18 (the CDIM, for example), is shown in the pseudo code of Figure 12. As described below, this routine is the starting point and the  
20 return point for many functions, described below with reference to Figures 13 - 22, that are performed by the system.

The first step, step 240, which is executed upon invoking this routine determines if one of the Incubators  
25 20 is out of its temperature specification or if another fault has placed the system into a non-functioning mode. If either of these conditions is detected, control is

- 35 -

transferred to the ALARM STATE routine, described below with reference to Figures 15a and 15b.

5 The next step, 242, examines the status of each of Incubators 20 to determine if any tray has completed its incubation period. If so, control is passed to the INCUBATION COMPLETE routine, described below with reference to Figures 14a and 14b.

10 If neither of the above conditions is detected, the operator is prompted, at step 244, to either enter a function from the menu or to wand the bar code identifier on a tray. Depending on the entered function, control is passed to the appropriate routine, described below. At step 246, the CONJUGATE RECONSTITUTION routine is described below with reference to Figures 21a and 21b. The BEAD  
15 ADDITION routine, invoked at step 248, is described below with reference to Figures 20a and 20b. At step 250, the OPD PREPARATION routine is described with reference to Figures 22a and 22b. The PRINT DI DATA routine is invoked at step 252 to print data provided by an Incubator 20.  
20 This routine is described with reference to Figure 16. The PRINT CDIM ERROR LOG routine is invoked at step 254 to print errors noted by the Incubator Monitor 18. This routine is described below with reference to Figure 17. At step 256, the COVER SEAL LOT NUMBER ENTRY routine,  
25 described below with reference to Figure 19, is invoked to add a new cover seal lot number to the system library. At step 258, the routine TIME OUT OF COOLER is invoked to allow a supply of conjugate beads to be removed from the

- 36 -

cooler for use by the system. This routine is described below with reference to Figure 18.

At step 260, no function was entered but a tray was wanded. In this instance, the MONITOR AND FUNCTION ENTRY routine invokes the BEGIN INCUBATION routine,  
5 described below with reference to Figures 13a and 13b.

The BEGIN INCUBATION routine is executed in Incubation Monitor 18. It starts at step 262, shown in Figure 13a. At this step, the technician is prompted to  
10 enter his identification code. At step 264, the system determines if the technician is qualified to initiate an incubation. If not, the network is notified at step 266 and control is returned to step 262 until an authorized identification code is received.

Otherwise, at step 268, the IM 18 asks network  
15 24 for instructions, or protocol requirements, on how this tray is to be handled. At step 270, the IM 18 verifies that the tray is ready to go into incubation based on the instructions received from the network. If IM 18  
20 determines, at step 272, that the tray is not ready, IM 18 displays an appropriate message and returns to the MONITOR AND FUNCTION ENTRY routine, described above with reference to Figure 12. If the tray is ready, then, at step 274, the IM 18 attempts to locate an available Incubator 20 slot.

The next step in the sequence is step 276 of  
25 Figure 13b. At this step, if a Incubator 20 exists which is in the proper mode, at the proper temperature and which

- 37 -

has an available slot, the IM 18 sets the timer for the slot to the time specified. If no Incubator 20 slot is available, the IM 18 gets the time for room temperature incubation from PC 16 by accessing the first row OPD  
5 dispense time stamp. In either instance, control is transferred to step 284.

At step 278, if an Incubator 20 is available but is not set up for the present protocol requirements for the tray, IM 18 will set the mode and temperature for the  
10 available Incubator 20 and will not allow Incubator 20 to be used until it has reached its set point. At step 280, IM 18 times-out the selected Incubator 20. If it does not reach its set point in a predetermined time limit, the ALARM STATE routine is invoked. Otherwise, when Incubator  
15 20 reaches its set point, IM 18 will set the timer for the slot and instruct the technician to insert the tray.

At step 284, if this is a conjugate or OPD incubation, IM 18 will prompt the technician to verify well volumes and place a cover seal on tray. At steps 286 and  
20 288, the technician either places the tray into the identified slot and pushes the button next to the slot or places it on the counter near a room temperature probe and pushes the ROOM TEMP button on IM 18 if the incubation is to be at room temperature. In response to this last step,  
25 IM 18 stores the incubation start time, Incubator 20 serial number, IM 18 software revision module the slot number, mode, temperature used and tray number in its internal memory. At step 292, IM 18 also passes this information to the DMS II through network 24. Finally, at step 294, IM

- 38 -

18 returns control to the MONITOR AND FUNCTION ENTRY routine, described above with reference to Figure 12.

5           Figures 14a and 14b describe the tray incubation complete routine. The first step in this routine, Step 298, IM 18 emits an alarm indicating that a tray needs to come out of its incubator. At step 300 the technician enters his identification code. At step 302, IM 18 verifies that the identification code is that of a qualified technician.

10           If, at step 304, the technician is not identified as qualified, IM 18 notifies the network of the unqualified technician status and returns control to step 298 until a qualified technician identification code is entered. If, at step 306, the technician is identified as  
15           being qualified, he or she is allowed to remove the tray from the Incubator 20 slot that is emitting the alarm. The next step in this process is step 308 of Figure 14b. At this step, the technician wands the bar code on the tray. At step 310, the technician pushes the button next to the  
20           slot where the tray was just removed or pushes the room temperature incubation button if the tray was being incubated at room temperature. At step 312 IM 18 receives the incubation stop time and from this calculates elapsed incubation time. At step 314 IM 18 stores the stop timer  
25           value for the tray number that was removed from the Incubator 20 slot. At step 316 IM 18 reports the stop time and elapsed time for the tray number to the network. At Step 318, IM 18 stores the minimum and maximum incubation temperatures and the mean temperature for the tray number.



- 39 -

At step 320 IM 18 notifies network 24 of these temperatures. Step 322 completes this routines by returning to the MONITOR AND THE FUNCTION routine.

5 The ALARM STATE routine is described with  
reference to figures 15a and 15b. At the first step in  
these figures, step 324, IM 18 goes into a severe alarm  
state which causes IM 18 to emit an alarm sound which is  
different than the alarm of Incubator 20. At step 326, IM  
18 stores an indication of the alarm, IM 18 serial number,  
10 and the serial number of Incubator 20 which caused the  
alarm.

Next, at step 328, IM 18 notifies network 24 of  
the type of alarm, either rotate, air flow, or temperature.  
IM 18 then sends its serial number and the serial number of  
15 Incubator 20 to network 24. At step 330, IM 18 determines  
whether another Incubator 20 is available to receive the  
trays in the failing Incubator 20. At step 332 if no other  
Incubator 20 is available then IM 18 voids all the  
associated trays and returns to the monitor and function  
20 entry routine, described above with reference to figure 12.  
Prior to returning to the monitor and function entry  
routine, IM 18 notifies network 24 of the voided trays at  
step 334.

At step 336, a substitute Incubator 20 has been  
25 found. Accordingly, IM 18 notifies the technician of which  
trays to move and where they are to be moved. At step 338,  
the technician is prompted to enter his identification  
code. At step 340 IM 18 verifies the identification code

- 40 -

against the network database qualified technicians. If, at step 342, the technician is not qualified to move the trays to the new incubator, IM 18 will notify the network of the unqualified technician status and transfer control to step 5 336 which will again prompt for a technician ID number until a qualified technician is found.

The ALARM STATE routine continues at step 344 of Figure 15b. At this step, once a qualified technician has been located, he is allowed to remove a specified tray from 10 the failing Incubator 20. At step 346, the technician wanders the tray to read its bar coded identification number. At step 348, the technician inserts the tray into the specified slot of the new Incubator 20 and wanders the tray position numbers. At step 350, the technician pushes the 15 tray slot button next to the specified slot. At step 352, IM 18 continues to monitor the tray in the new documented position and stores an additional IM 18 serial number, Incubator 20 serial number IM 18, software revision number, slot number, mode, temperature and tray number in its local 20 memory. At step 354 IM 18 notifies network 24 of this data. At step 356 control is returned to step 344 which repeats the above steps until all the trays have been moved from the failing Incubator 20 to the newly found Incubator 20. At step 358 IM 18 returns to the function which 25 invoked the alarm state routine or returns to the MONITOR AND FUNCTION ENTRY routine.

Figure 16 describes the function performed by the PRINT INCUBATOR 18 DATA routine. As described above, this routine prints data stored in a specified Incubator

- 41 -

20. The first step in this routine, step 360, is to prompt  
the technician to specify which Incubator 20 the data is to  
be read from. At step 362, the technician enters this  
Incubator 20 identification number. At step 364, IM 18  
5 prints to a printer, temperature, mode, and assay by tray  
slot and the time that the tray is due to come out for each  
tray in the dynamic integrator. At step 366, IM 18 returns  
to the MONITOR AND FUNCTION ENTRY routine. Figure 17  
describes the PRINT CDIM ERROR LOG routine. At step 368,  
10 this routine causes IM 18 to print the error log which  
includes at least the last 100 error codes, to a printer.  
At step 370 IM 18 returns control to the FUNCTION AND DATA  
ENTRY routine described above with reference to figure 12.

Figure 18 describes the time out of cooler  
15 routine. In the first step of this routine, 372, IM 18  
prompts the technician to enter his or her identification  
code. At step 374, IM 18 verifies the entered  
identification code against the database of qualified  
20 technicians. At step 376 if IM 18 determines that the  
technician is not qualified, it notifies network 24 that  
an unqualified technician has attempted access and returns  
to step 372 to request the technician to enter a valid ID.

Once a valid technician ID has been entered, at  
step 378, the technician is prompted to wand the lot number  
25 of the conjugate or bead container which is to be removed  
from the cooler. At step 380, IM 18 records the time at  
which the conjugate bead container was removed from the  
cooler and stores this and the associated lot number in its  
local memory. At step 382, IM 18 notifies network 24 of

- 42 -

this data. Next, at step 384, IM 18 polls network 24 to determine the validity of this lot number. At step 386 if the lot number is determined to be invalid, then IM 18 displays a message to that effect and returns control to  
5 step 378 to prompt the technician to wand another container. At step 388, IM 18 sets a timer on network 24 which prohibits the use of the removed container until 30 minutes after the container has been removed to allow the conjugate or beads to come up to room temperature. At step  
10 390, the routine returns control to the MONITOR AND FUNCTION ENTRY routine.

Figure 19 shows the program flow for the COVER SEAL AND LOT NUMBER ENTRY routine. The first step in this routine, step 392, prompts the technician to enter his or  
15 her identification code. At step 394, IM 18 verifies the entered code against the database of qualified technicians. If, at step 396, IM 18 determines that the technician is not qualified it notifies the network of the unqualified technician status and returns to step 392 which repeats the  
20 above process until a valid technician ID is entered. At step 398, once a valid technician has been located, the technician is prompted to enter the cover seal lot number. At step 400, IM 18 polls network 24 to determine the validity of this lot number. If, at step 402, the lot  
25 number is invalid IM 18 displays a message to that effect and returns control to step 398 to request the technician to enter another cover seal and lot number. If at step 403, the lot number is found to be valid, then IM 18 causes the lot number to be entered into the network library for  
30 use by system 10. At step 404, the COVER SEAL LOT NUMBER

- 43 -

ENTRY routine returns control to the MONITOR AND FUNCTION ENTRY routine, described above with reference to figure 12.

Figure 20a describes the program flow for the BEAD ADDITION routine. At the first step of this routine, step 406, the technician is prompted to enter his identification code. At step 408 IM 18 verifies that the Identification code is valid against the network database of qualified technicians. At step 410, if it was determined that the technician is not qualified, IM 18 notifies network 24 of the unqualified technician status and returns to step 406 until a valid technician identification code is entered.

At step 412, once a valid technician identification code has been entered, the technician is prompted to wand the Identification code number into IM 18 from the tray. At step 414, IM 18 requests from network 24, the processing protocol to be used for this tray. At step 416, the IM 18 verifies that the tray is ready for bead addition based on the received protocol sequence. Next, at step 418, if the IM 18 determines that the tray is out of sequence, and therefore not ready for bead addition, it notifies the technician as to where the tray should be in the system and transfers control to the MONITOR AND FUNCTION ENTRY routine.

Otherwise, if the tray is ready for bead addition, the technician is prompted to enter the bead lot number at step 420. After the technician has entered the lot number, IM 18 polls network 24, at step 422, to

- 44 -

determine the validity of the lot number. If IM 18 determines that the lot number is invalid, then, at step 424, it displays a message to that effect and returns control to step 420 until a valid lot number is entered.

5 At step 426, once a valid lot number has been entered, it is recorded and displayed by IM 18.

At step 436, IM 18 displays a message requesting the operator to dispense the beads, verify that the beads

10 have been dispensed, tap the tray, cover the tray and press "done" when the tray has been covered. Next, at step 437, the technician (operator) is prompted to verify the lot number against the display which shows the current library lot number. If the lot number of the cover seal is not the

15 one displayed, the technician may enter a new lot number into the library. At step 438 in the process, the technician dispenses the beads into the tray, visually inspects the tray and taps the tray to settle the beads. At step 440, the technician places the cover seal on the

20 tray. At step 442, when the technician has pressed the "done" key on IM 18, IM 18 stores the time the beads were dropped and stores the cover seal lot number information. At step 444, IM 18 notifies network 24 of the time the beads were dropped and the cover seal lot number. At step

25 446, the BEAD ADDITION routine transfers control the BEGIN INCUBATION routine, described above with reference to figures 13a and 13b.

- 45 -

Figure 21a and 21b show the program flows for the CONJUGATE RECONSTITUTION routine. This routine starts at step 448 when IM 18 prompts the technician to enter his identification code. At step 450, IM 18 verifies the technician identification code against network 24 database of qualified technicians. If, at step 452, IM 18 determines its technician is not qualified, it notifies network 24 of the unqualified technician status and returns control to step 448 until a valid technician identification code is entered.

At step 454, the technician is prompted by IM 18 to wand the identification information from the conjugate concentrate and conjugate diluent containers. At step 456, IM 18 receives this identification information and polls network 24 for its validity. Network 24 responds with the validity data at step 458. If, at step 460, IM 18 was informed by network 24 that the conjugate diluent or conjugate concentrate were invalid, it displays an error message and returns to step 454 until information from valid containers is entered. At step 462, when valid containers are available, the technician is prompted to combine the concentrate and conjugate diluent.

Next, at step 464, the technician defaces and relabels the bottle containing the diluted concentrate. At step 466, a new conjugate bar code is placed on the bottle containing the diluted concentrate. At step 468, the technician is prompted to wand the new conjugate bar code. IM 18 receives this information and passes it on to network 24.

- 46 -

Network 24 responds at step 470 by providing a time and date at which the conjugate was reconstituted and an expiration time and date. At step 472, IM 18 stores the conjugate lot number and expiration time and date, the  
5 diluent lot number and expiration time and date, the time and date at which the reconstitution was started, and the lot number and expiration time and date for the diluted conjugate. At step 474, IM 18 passes this information on to network 24.

10 Next, at step 476, the technician is prompted to add the displayed information to the label on the bottle that holds the reconstituted conjugate. At step 478, IM 18 starts a 30-minute timer. This timer defines a time in  
15 which the conjugate is allowed to stabilize before it may be used. At step 480, when the 30-minute timer expires, IM 18 stores the time and date at which the reconstitution was completed. At step 482, IM 18 notifies network 24 of this time and date. In the exemplary system, network 24 will  
20 not allow the diluted conjugate to be used before the 30-minute timer has elapsed or after its expiration time and date. At step 486, the final step of the CONJUGATE RECONSTITUTION routine, control is transferred back to the MONITOR AND FUNCTION ENTRY routine.

25 Figures 22a and 22b describe the program flow for the OPD PREPARATION routine. This routine is very similar in its operation to the CONJUGATE RECONSTITUTION routine, described above with reference to Figures 21a and 21b.



- 47 -

At the first step of the OPD PREPARATION routine, step 488, IM 18 prompts the technician to enter his identification code. At step 490, IM 18 verifies the entered identification code against network 24 database of qualified technicians. If, at step 492, IM 18 determines that the entered technician identification code is not valid or that the technician is not qualified, it notifies network 24 of unqualified technician status and returns control to step 488 until the identification code of a qualified technician is entered.

At step 494, IM 18 prompts the technician to wand the identification information from the OPD tablet container and the diluent container. At step 496, IM 18 passes this information on to network 24 to determine the validity of the tablets and diluent. At step 498, network 24 responds with this validity information. If, at step 500, IM 18 is informed that the OPD tablets or the diluent is invalid, it displays an error message and returns control to step 494 until information from valid containers is entered.

At step 502, the OPD tablets and diluent are valid and the technician is prompted to combine the tablets with the diluent. At step 504, the technician is prompted to deface and relabel the bottle which holds the OPD substrate (which is a combination of tablets and diluent), and at step 506, to add a new OPD substrate bar code. At step 508, the technician wands the new substrate bar code into IM 18 which passes it on to network 24. Next, at step 510, network 24 responds by sending the date and time at

- 48 -

which the OPD was prepared and an expiration date and time for the OPD solution to IM 18.

5 The next step in the process is step 512 of Figure 22b. At this step, IM 18 stores the OPD lot number, expiration time and date, diluent lot number with its expiration time and date, the time and date at which the reconstitution was started, lot number of the prepared OPD substrate and its expiration date. At step 514, IM 18 passes this information on to network 24.

10 At step 516, the technician is prompted to add this information to the OPD substrate label. Next, at step 518, IM 18 starts a 30-minute timer to allow the newly prepared OPD solution to stabilize. At step 520, when the timer expires, IM 18 stores the time and date at which the reconstitution was completed. Next, at step 522, IM 18  
15 notifies network 24 of the time and date at which the OPD reconstitution was completed. At step 526, the last step of the OPD PREPARATION routine, control is returned to the monitor and function entry routine described above with  
20 reference to Figures 12a and 12b. System 10, through network 24 database, will not allow the OPD substrate to be used during the 30-minute stabilization period or after its expiration time and date.

25 Figure 23 is a flow chart diagram which illustrates an exemplary test which may be performed using the system shown in Figure 1.

- 49 -

In the first step of this test, step 570, the reagents are prepared. In the exemplary test three reagents are used, a conjugate, OPD and an acid. Depending on the status of system 10, these reagents may be already available at room temperature in which case step 570 simply consists of attaching an external pump device to each of the bottles. If the reagents are prepared but stored in the cooler, they must be removed and allowed to stand at room temperature for 30 minutes or some other appropriate time, before they may be used. This process is described above with reference to Figure 18.

If no reagents are available at room temperature and none are in a prepared state and stored in the cooler, then reagents must be prepared; conjugate from concentrate and OPD from OPD tablets. The processes for preparing each of these reagents are described above with reference to Figures 21A and 21B for the conjugate, and with reference to Figures 22A and 22B for the OPD.

The end result of step 570, which is described below in greater detail with reference to FIG. 24, is a supply of three reagents, a conjugate, an OPD solution and an acid.

With the reagents available and at room temperature, the next step in the process, step 572, is to process a tray of samples through a processing unit 16 using the conjugate. The steps for processing each tray using a reagents are described below with reference to Figures 25, 26 and 27A and 27B.

- 50 -

Once the tray samples has been processed using the conjugate at step 572, the next step, 574 processes the tray through bead addition and incubation stages. In the exemplary embodiment of the invention both of these  
5 processing stages are performed on Incubator Monitor 18, and are described above with reference to Figures 20A and 20B for the bead addition and Figures 13A, 13B, 14A and 14B for the incubation.

The next step in the process shown in Figure 23,  
10 step 576, processes the tray of samples using the OPD solution as the reagent. As with step 572, the steps performed by Processing Center 16 to implement step 576, are described below with reference to Figures 25, 26, 27A and 27B.

15 After the tray has been processed to Processing Center 16 using the OPD, the next step in the process is an OPD incubation performed at step 578. The incubation step is performed at Incubator Monitor 18, and Incubators 20 or 21 as described above with reference to Figures 13A, 13B,  
20 14A and 14B.

When the tray has been removed from the incubator after its OPD incubation, it is processed again though Processing Center 16 using the acid as the reagent. This is illustrated by step 580 and described below with  
25 reference to Figures 25, 26, 27A and 27B.

At the end of processing the tray, using the acid as the reagent, the processing center 16 reads the

- 51 -

results of the test at step 582. The manner in which the test described above with reference to Figure 23 is implemented is described below in greater detail with reference to Figure 31.

5                   Figure 24 is a program flow diagram which illustrates the PREPARE REAGENTS routine shown in Figure 23. This routine assumes that prepared reagents are available in a cooler (now shown). The next step in this routine, step 586 is to bring the reagents to the work station and, at step 588, to allow the reagents to come to room temperature for 30 minutes. In the exemplary embodiment of the invention these two steps are performed using the incubator monitor 18; and are described above with reference to Figure 18.

10                   The next step in the process shown in Figure 24, is to put the TRI-CONT onto each of reagent bottles. In the exemplary embodiment of the invention, the TRI-CONT is a small external pump mechanism which is attached to each of the reagent bottles, so that, when the bottles are fitted into a Processing Center 16, Processing Center 16 will be able to dispense the reagent from the bottles. Step 592 completes the prepare reagents routine by priming each of the TRI-CONTS to make them ready for use by their Processing Centers 16.

15                   Figure 25 is a program flow diagram which illustrates the operation of a valid tray routine which is one of the routines run to process a tray through one of the Processing Centers 16. At the first step of this

- 52 -

routine. step 594, Processing Center 16 prompts the operator to put the tray into Processing Center 16 entrance chute. At step 596, Processing Center reads the tray bar code or requests keyboard entry of the bar code by the technician. Next, at step 598, Processing Center 16 determines if the tray is active and has valid control samples. At step 600, if Processing Center 16 receives information from network 24 indicating that the tray is invalid, it displays an appropriate message and instructions to the operator and returns to step 594 to wait until a valid tray is entered.

Figure 26 is a program flow diagram which illustrates the operation of the VALID TECHNICIAN routine. This is a second routine run in processing the trays as in steps 572, 576 and 580 of Figure 23. At the first step of this routine, step 602, Processing Center 16 displays a message asking the technician to enter his ID code. At step 603 the operator enters his ID code and at step 604 Processing Center 16 polls the system by network 24 to determine if the technician ID is valid. At step 605, if the technician ID is found not to be valid the processing center displays an appropriate error message and instructions and returns control to step 602 until a proper technician ID is entered.

Figures 27A and 27B are flow diagrams which illustrate the bulk of the processing performed by Processing Center 16 in processing a tray with one of the reagents (i.e. the conjugate solution, the OPD solution or the acid). At the first step in this process, Processing

- 53 -

Center 16 displays a message asking the operator to enter the number of trays in the batch. At step 608, the operator enters the number of trays. At step 610, Processing Center 16 displays the assay name and number  
5 reagent location and asks the operator to enter the reagent identification information from the bar code on the reagent bottle. At step 612, the operator enters the reagent information using, for example, bar code reader 19. Alternatively, the operator may use a dedicated bar code  
10 reader (not shown), coupled to Processing Center 16.

Upon receiving the reagent identification information, at step 614, Processing Center 16 polls system 10 through network 24 to determine if the reagent is valid. At step 616, Processing Center receives a message from  
15 network 24 indicating the status of the reagent.

If, at step 618, the reagent is invalid, Processing Center 16 displays an appropriate message and loops back to step 610 until a valid reagent lot number is entered.

20 At the next step of the process, step 620, Processing Center 16 displays a message asking the operator to enter the volume of the reagent. At 622, the operator reads the approximate volume from the reagent bottle using demarcations on the bottle. At step 624 the operator  
25 enters this volume via the keyboard (now shown) coupled to Processing Center 16. Next, at 626, Processing Center 16 displays a number of samples and controls which can be processed and requests entry ("yes" or "no") if the number

- 54 -

of samples and controls in the batch is less than or equal to the number that is displayed. At step 628, if the answer entered by the operator is "no", Processing Center 16 displays an appropriate message and loops back to step 5 620 until the operator enters an volume which is suitable for the number of samples.

The next step in the process is step 630 of Figure 27B. At this step Processing Center 16 displays a message asking the operator to load the reagent on a 10 particular station number and enter the reagent station number. At step 632, the operator loads the reagent into the displayed station including tubing and tip necessary to dispense the reagent. At step 634, the operator enters the station number via the bar code reader 19 or the bar code 15 reader (not shown) next to the station, or alternatively by the keyboard (not shown) coupled to Processing Center 16. Next, at step 636, Processing Center 16 displays the message "valid station." At step 638, the operator enters "yes" or "no" to verify that the station is valid. At step 20 640, if the operator answered "no" to the valid station request, Processing Center 16 displays an appropriate message and loops back to step 630 until the reagent container is loaded into a valid station. At step 642, Processing Center 16 displays the message "insert tray." 25 At step 644, Processing Center 16 displays a message "ready." At step 646, the operator enters either "yes" or "no" to approve processing. If, at step 648, the operator has entered a "no", Processing Center 16 displays an appropriate message to remove the tray.



- 55 -

At step 650, once the reagent is loaded, Processing Center 16 pulls the tray in to dispense the reagent. At step 652, Processing Center 16 saves the time when the processing begins. At step 654 Processing Center  
5 16 completes processing and at step 656, passes the tray to the exit chute after the reagent has been dispensed. Next, at step 658, Processing Center 16 sends a processing complete message, the dispense start time and stop time, and any error conditions which were encountered and any  
10 voided wells to the Data Management System 12 via network 24. At step 660, Processing Center 16 emits a beep to signal the operator that it is time to remove the tray from the exit chute if necessary. At step 662, the operator removes the tray from the exit chute.

15 When the assay requirements as defined by the assay protocol or "recipe" have been completed, the results will be collated by DMS 12. Figure 28, step 670 describes accepting these results. The technician may choose to either accept the batch results or not. If the results are  
20 accepted, then at step 672 the results go through the decision making tables stored in the database of DMS 12. At step 674, the technician may then designate whether the results are to be used in statistical reports. This allows the user to screen invalid or erroneous runs out of reports  
25 used for analysis.

Subsequent to running various tests on various samples, the user has options to generate reports. Figure 29 describes some of these options. Operational reports include batch reports on current sample results such as

- 56 -

those described at step 676. Step 678 describes an exception report, that is, a list of all unacceptable sample results. Step 680 describes the retest report which details those samples requiring further iterations of tests based on statistical accuracy measurements. Step 682 describes typical Quarantine Reports which will provide the user a list of sample identification numbers for those samples reactive to a particular test.

The user may use this report to quarantine materials associated with the biological material being tested. Management or decision support reports typically include comprehensive batch records and statistical reports such as described in figure 30. Step 684 describes a batch record which will include a history of all results, decisions, processing steps, reagents and commodities used (i.e. lot number, expiration date), technicians performing steps, instruments used for processing, times processing took place, and other parameters. Step 686 describes statistical reports which will be defined by the user. Examples may include reports describing control performance, reagent performance, incidence of reactive results and the like.

All of this information (data) can be stored in a long-term storage device (not shown) on or off network for later retrieval.

There is shown in Figure 23 a flow diagram for the operation of the method of the present invention. There is shown in block 550 the step of identifying the

- 57 -

sample to be tested. In the exemplary embodiment, this involves scanning the bar code on the test tube containing the sample to be tested. In block 552, the identification information of the sample read in block 550 is communicated  
5 to the storage device and database which contains the information regarding the test(s) to be undertaken for the particular sample. This information will be communicated back to Pipettor 14 or FPC, along network 24 to pipet the proper aliquot to the proper test sites.

10 In block 54, the particular test site is identified. The identification is completed by scanning the bar code on the test tray containing the test sites. The test sites are set up as a matrix with their own  
15 particular identification and database coordinated to the bar code label identifying the test tray. In block 556 the identification information for the test site (test tray) is communicated over network 24 to the storage device and  
20 database containing this information. It is possible that this information is stored locally at a pipettor 14 or at a pipettor monitor 13. The feature of the storage is that the information is accessible to the pipettor and any other device which may require access to this information.

In block 558, the aliquot of sample is delivered to the appropriate test site, depending upon the test  
25 required for the sample. This information is provided by DMS 12 along network 24 to pipettor 14 (or pipettor monitor 13).

- 58 -

In block 560, the aliquot of sample is tested. This testing may involve reagent addition through a processing center 16 and/or incubation in an incubator 20. Further, testing may involve the addition of a coated bead or similar reagent delivery method. Further, testing may involve multiple passes through processing center 16 and/or incubation in an incubator 20. The exact test steps depend upon the test being conducted on the particular aliquot of sample. In block 562, results of the tests are read. In the exemplary embodiment, 16 is capable of reading the test result(s) by spectrophotometric variation of the test sample during the test process. Finally, in block 564 the test results are communicated to a memory storage device such as DMS 12 for later retrieval by the user of the system.

The present invention has been described in terms of exemplary embodiments. It is contemplated, however, that it may be practiced with modifications, some of which are outlined above, within the scope of the appended claims.

The Invention Claimed is:

- 1                   1. A system suitable for use in testing  
2 biological fluids, comprising:
- 3                   a data communications network;
- 4                   means for receiving a sample of a biological  
5 fluid having an identifying label and for communicating  
6 identifying information contained on the label to the  
7 network;
- 8                   means, coupled to the data communications  
9 network, for receiving the identifying information to  
10 determine at least one test to be performed on the sample;
- 11                   means, coupled to the data communications  
12 network, for testing the sample according to the determined  
13 test and for communicating results of the test to the  
14 network; and
- 15                   means for receiving and recording the results of  
16 the test from the network.
- 1                   2. A system for testing a sample of a biological  
2 fluid comprising:
- 3                   pipetting means for transferring a portion of  
4 the sample to a test site;

- 60 -

5 reagent testing means for applying at least one  
6 reagent to the test site;

7 means for incubating the test site, treated by  
8 the reagent testing means;

9 a network which is coupled to provide data  
10 communications among the pipetting means, the reagent  
11 testing means and the means for incubating; and

12 a data processing system, coupled to the  
13 network, for coordinating and tracking operations performed  
14 in the testing system by providing commands to and  
15 receiving data from each of the pipetting means, the  
16 reagent testing means and the means for incubating, via the  
17 network.

1 3. A method for testing biological fluids  
2 comprising the steps of:

3 identifying a sample of a biological fluid for  
4 testing;

5 communicating identification information of said  
6 sample to a first memory storage location on said  
7 communication network;

8 identifying a test site;

- 61 -

9                   communicating identification information of said  
10 test site to a second memory storage location on said  
11 communication network;

12                   delivering a predetermined amount of said sample  
13 to a test site;

14                   testing said sample;

15                   reading a result of said test;

16                   communicating said result to a third memory  
17 storage device, said third memory storage location on said  
18 communication network.

1                   4. The method of claim 3 wherein said step of  
2 testing said sample further comprises the steps of:

3                   transporting said test site to a test device;

4                   identifying said test site at said test device;

5                   communicating identification information of said  
6 test site at said test device to a fourth memory storage  
7 device, said fourth memory storage location on said  
8 communication network;

9                   communicating testing instructions to said test  
10 device over said communication network;

- 62 -

11                   adding reagents to said test site, according to  
12                   said testing instructions.

1                   5. The method of claim 4 further comprising the  
2                   steps of:

3                   transporting said test site to an incubator;

4                   identifying said test site at said incubator;

5                   communicating identification information of said  
6                   test site at said incubator to a fifth memory storage  
7                   location, on said communication network;

8                   communicating incubation instructions for said  
9                   sample to said incubator over said communication network;

10                  incubating said sample in said incubator  
11                  according to said incubation instructions.

1                   6. The method of claim 5 further comprising the  
2                   steps of:

3                   entering reagent information of reagents into a  
4                   memory storage device;

5                   identifying said reagents to be loaded onto said  
6                   test device;

7                   loading said reagents onto said test device;



8 reading said reagent information prior to  
9 conducting a test; and

10 determining whether said reagents are adequate  
11 for conducting said test.

1 7. A method for testing biological fluids  
2 comprising the steps of:

3 identifying a sample of a biological fluid for  
4 testing;

5 communicating identification information of said  
6 sample to a memory storage location, on a communication  
7 network;

8 identifying a test site;

9 communicating identification information of said  
10 test site to said memory storage location;

11 communicating delivery information for said  
12 sample to a delivery device over said communication  
13 network;

14 delivering said sample to said test site  
15 according to said delivery information;

16 transporting said test site to a reagent testing  
17 device;

18                    identifying said test site at said reagent  
19    testing device;

20                    communicating identification information of said  
21    test site at said testing device to said memory storage  
22    location;

23                    communicating test information to said testing  
24    device over said communication network;

25                    applying reagents to said test site according to  
26    said testing information;

27                    transporting said test site to an incubation  
28    device;

29                    identifying said test site at said incubation  
30    device;

31                    communicating identification information of said  
32    test site at said incubation device to said memory storage  
33    location; and

34                    communicating incubation initiation instructions  
35    to said incubation device over said communication network.

11/25

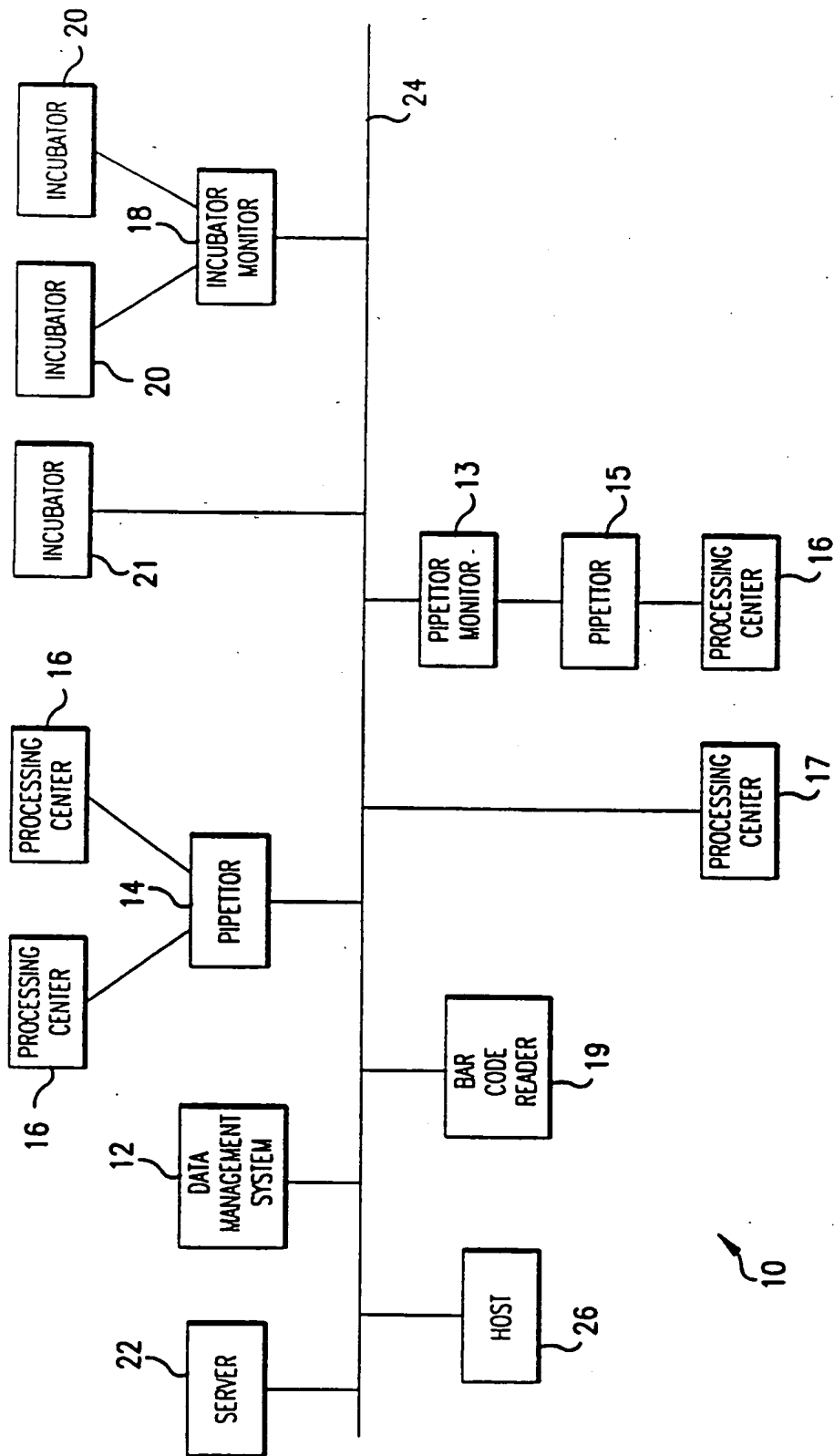


FIG.1

2/25

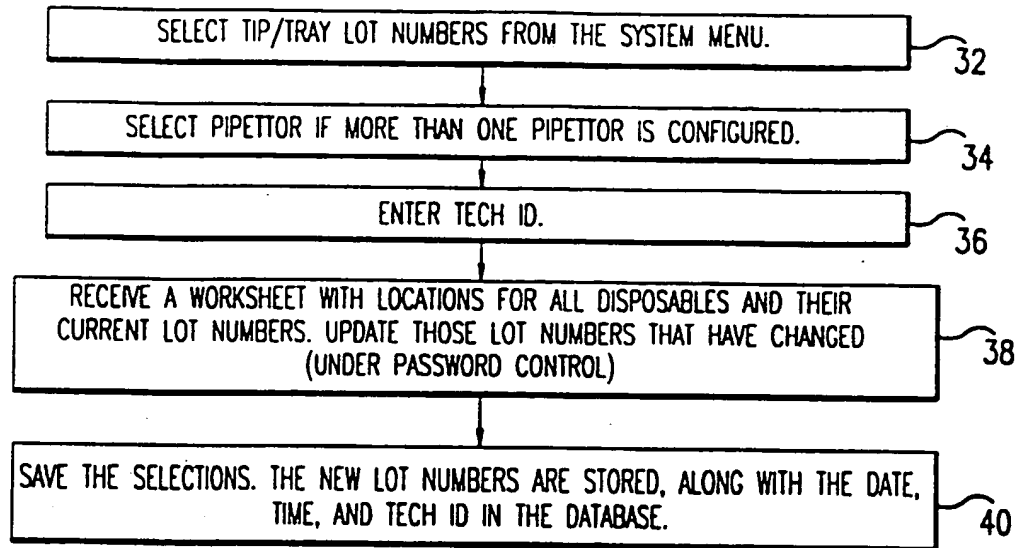


FIG. 2

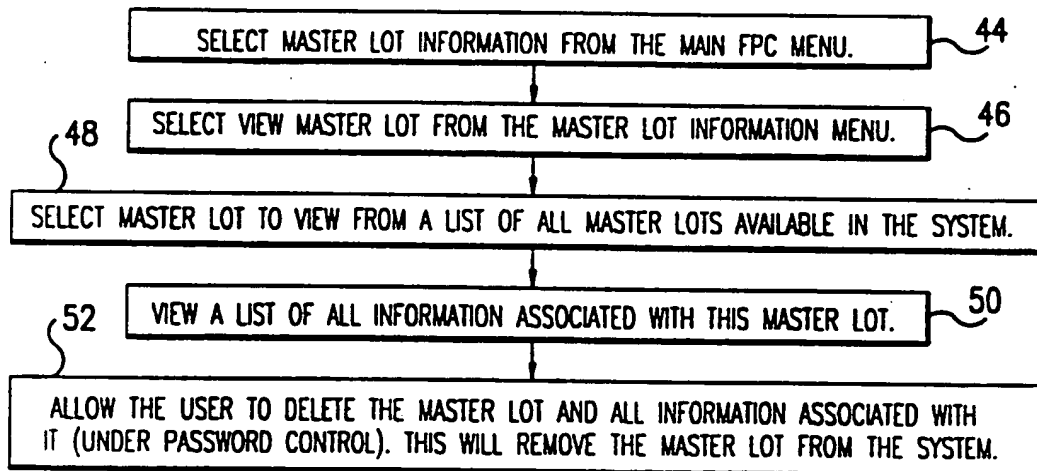


FIG. 3

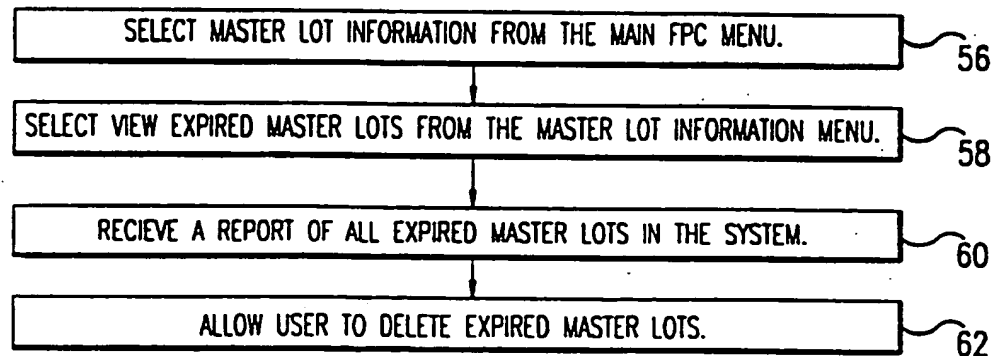


FIG. 4

3/25

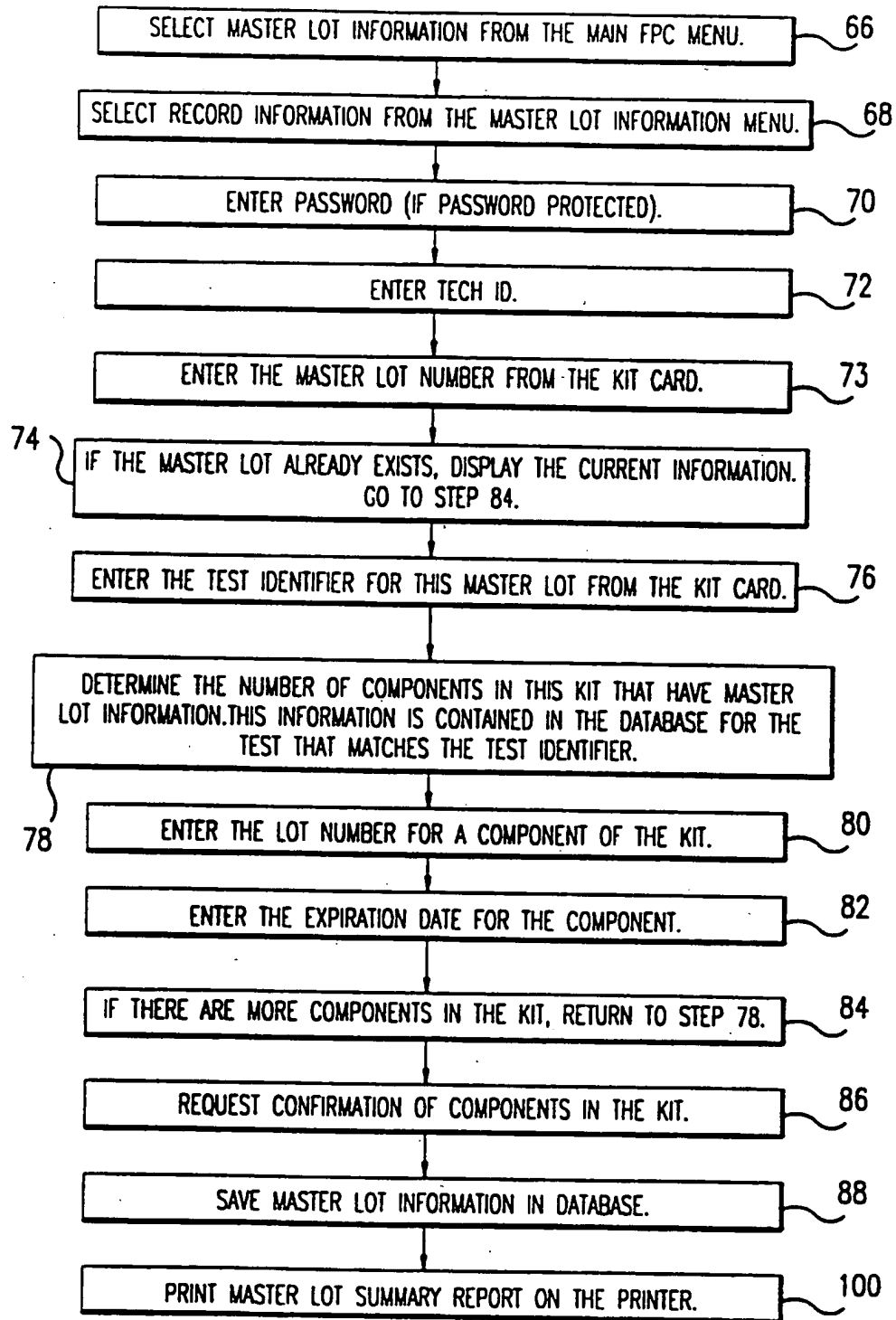


FIG.5

4/25

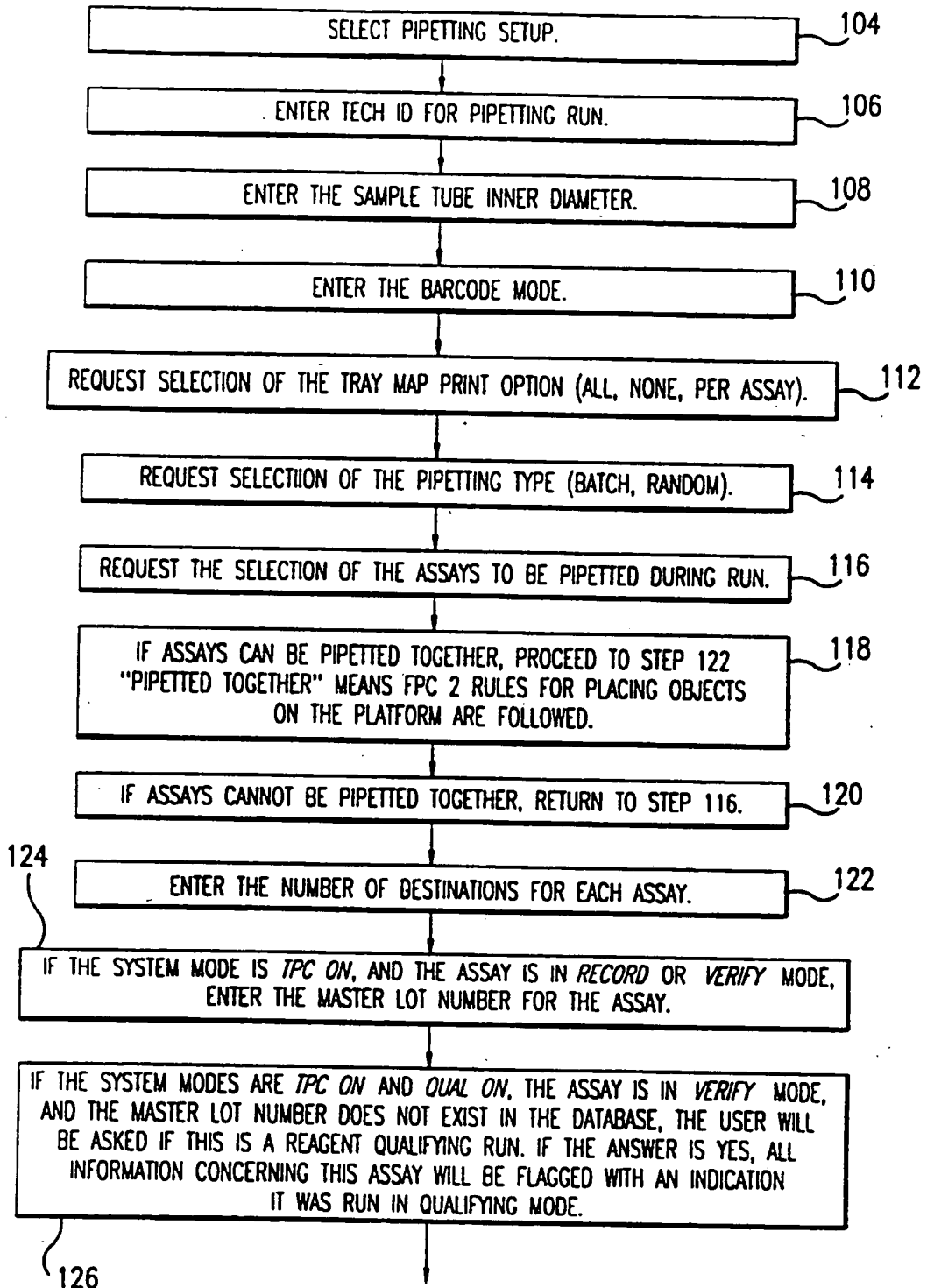


FIG. 6

5 / 25

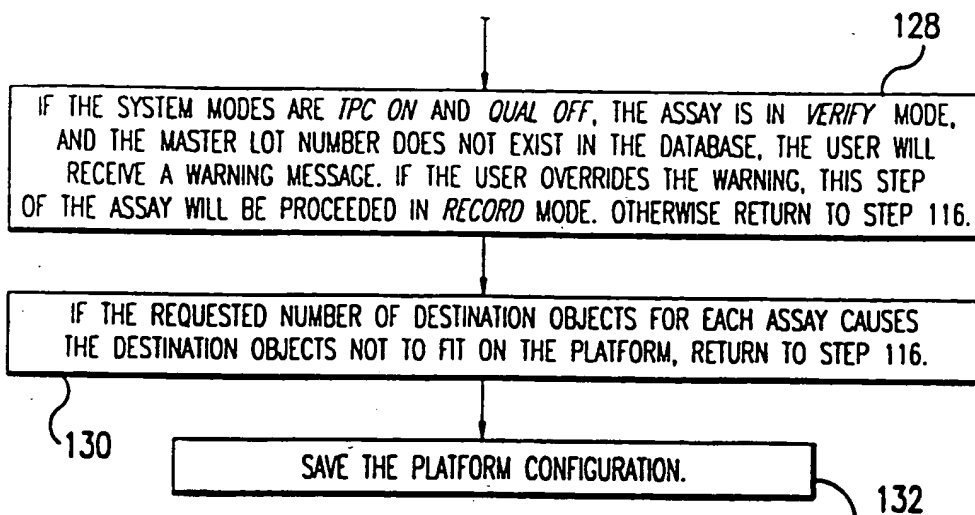


FIG.6 Cont.

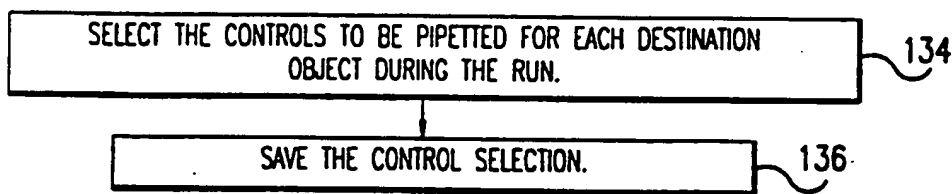


FIG.7

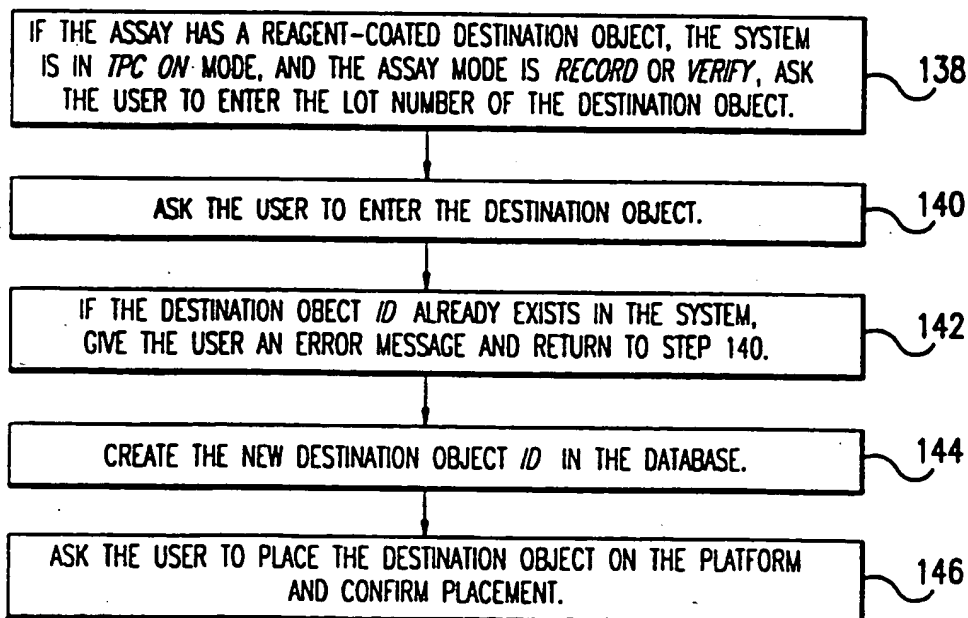


FIG.8

6/25

IF THERE ARE MORE DESTINATION OBJECTS THAT NEED TO BE PLACED, RETURN TO STEP 138. OTHERWISE, CONTINUE WITH STEP 138. 148

FIG.8 Cont.

DISPLAY A SCREEN INDICATING THE LOCATION FOR ALL REAGENTS AND CONTROLS ON THE PLATFORM. 150

IF THE SYSTEM MODE IS *TPC ON*, AND THE ASSAY MODE IS *RECORD* OR *VERIFY* MODE, REQUEST A LOT NUMBER FOR THE DILUENT OR CONTROL BOTTLE. TO OVERRIDE *RECORD* MODE, THE USER PRESSES <ENTER>. TO OVERRIDE *VERIFY* MODE, THE USER MUST CONFIRM A WARNING MESSAGE.

152 CHECK THE *ID* OF THE CONTROL AGAINST THE MASTER LOT LIST (*TPC ON* AND *VERIFY* MODE ONLY). 154

IF THE *ID* DOES NOT MATCH AN ENTRY IN THE LIST, GIVE THE USER AN ERROR MESSAGE AND RETURN TO STEP 152 (*TPC ON* AND *VERIFY* MODE ONLY) 156

IF THE EXPIRATION DATE OF THE CONTROL HAS BEEN REACHED, GIVE THE USER AN ERROR MESSAGE AND GO TO STEP 168 (*TPC ON* AND *VERIFY* MODE ONLY) 158

IF THE MASTER LOT CONTROL IS NOT CONSISTENT WITH THE OTHER CONTROLS IN USE FOR THIS ASSAY, GIVE THE USER AN ERROR MESSAGE AND GO TO STEP 168 (*TPC ON* AND *VERIFY* MODE ONLY) 160

ASK THE USER TO PLACE THE CONTROL BOTTLE IN THE SELECTED POSITION. 162

ASK THE USER TO BARCODE THE POSITION OF THE CONTROL BOTTLE. IF THE POSITION DIFFERS FROM THE EXPECTED POSITION, GIVE THE USER AN ERROR MESSAGE AND RETURN TO THIS STEP. 164

CONTINUE WITH STEP 152 FOR ALL REMAINING DILUENTS AND CONTROLS FOR THIS PIPETTOR RUN. IF ALL DILUENTS AND CONTROLS HAVE BEEN ASSIGNED POSITIONS, GO TO STEP 174. 166

FIG.9



7/25

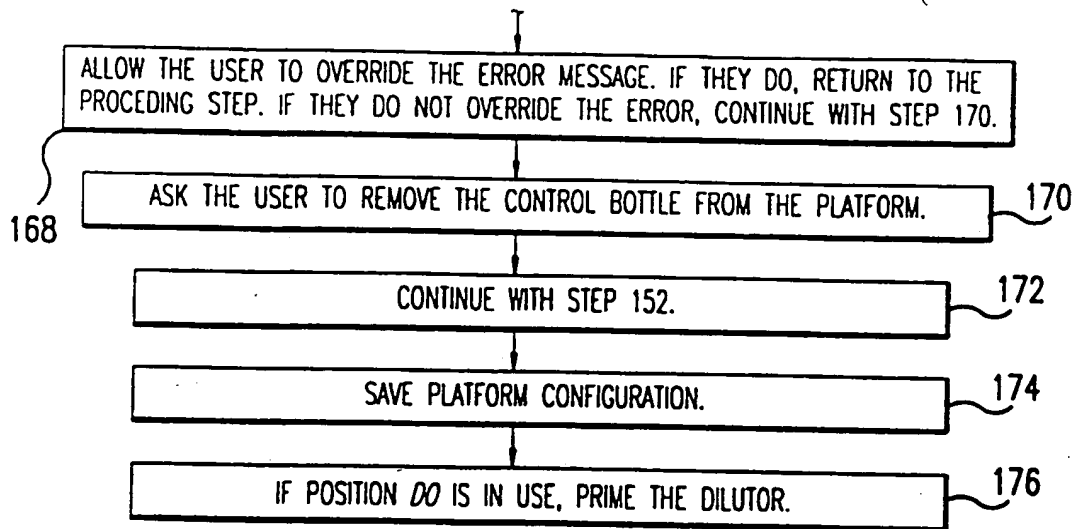


FIG.9 Cont.

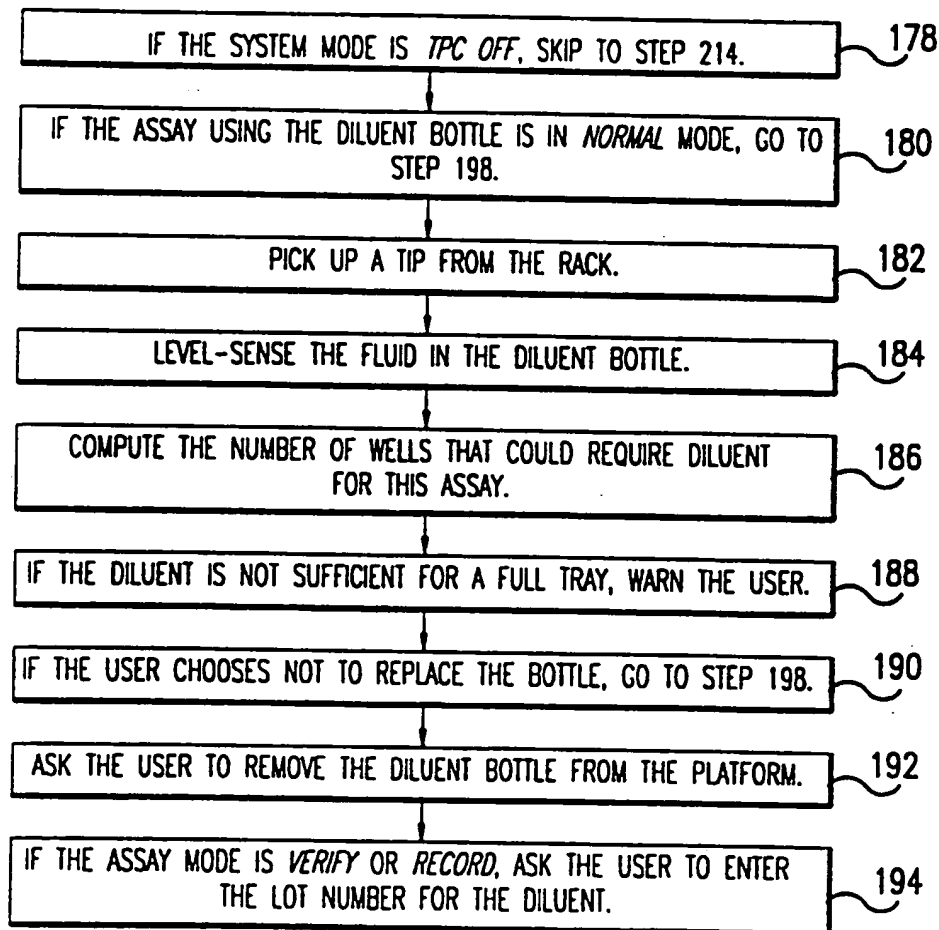


FIG.10

8/25

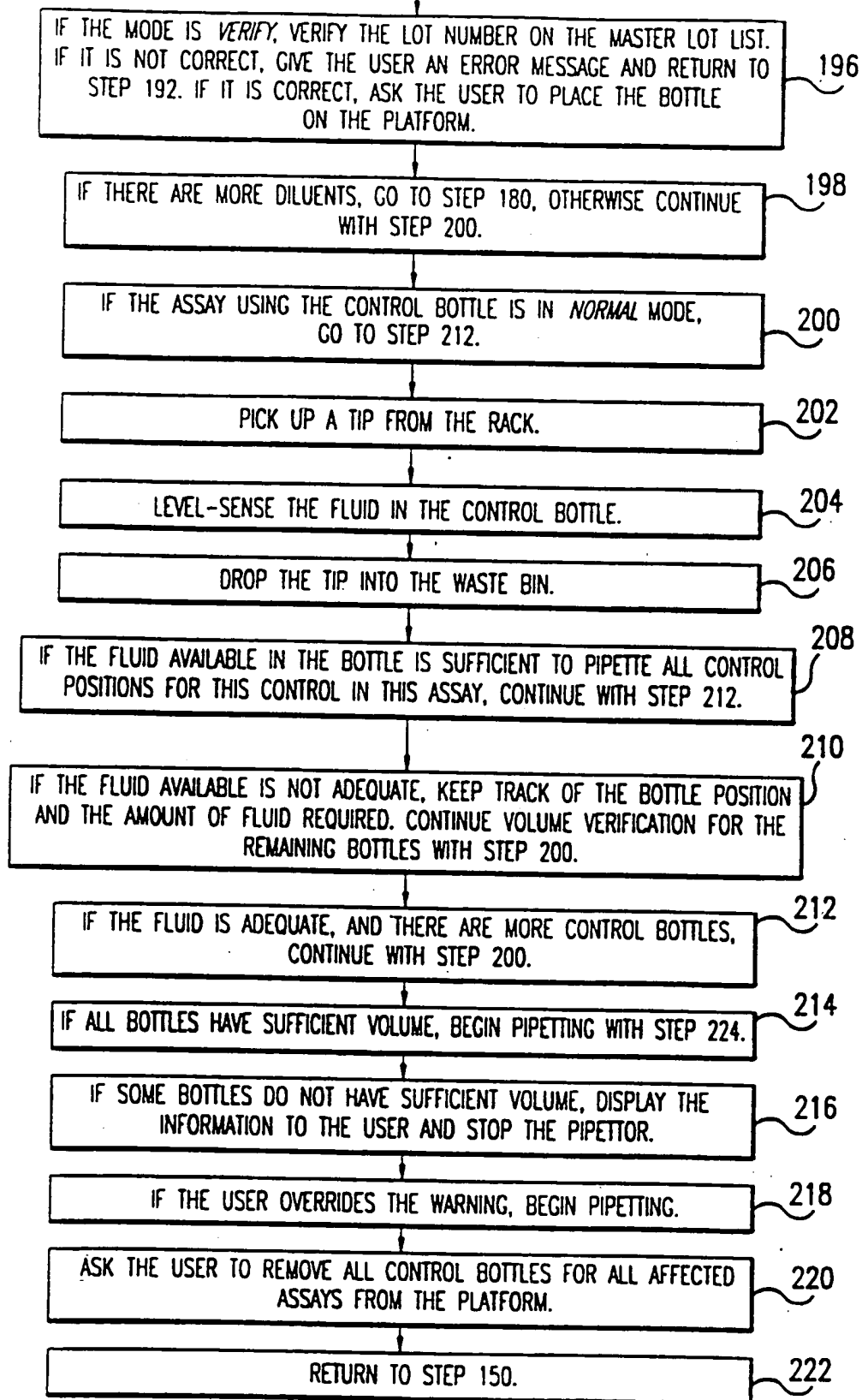


FIG.10 Cont.

9/25

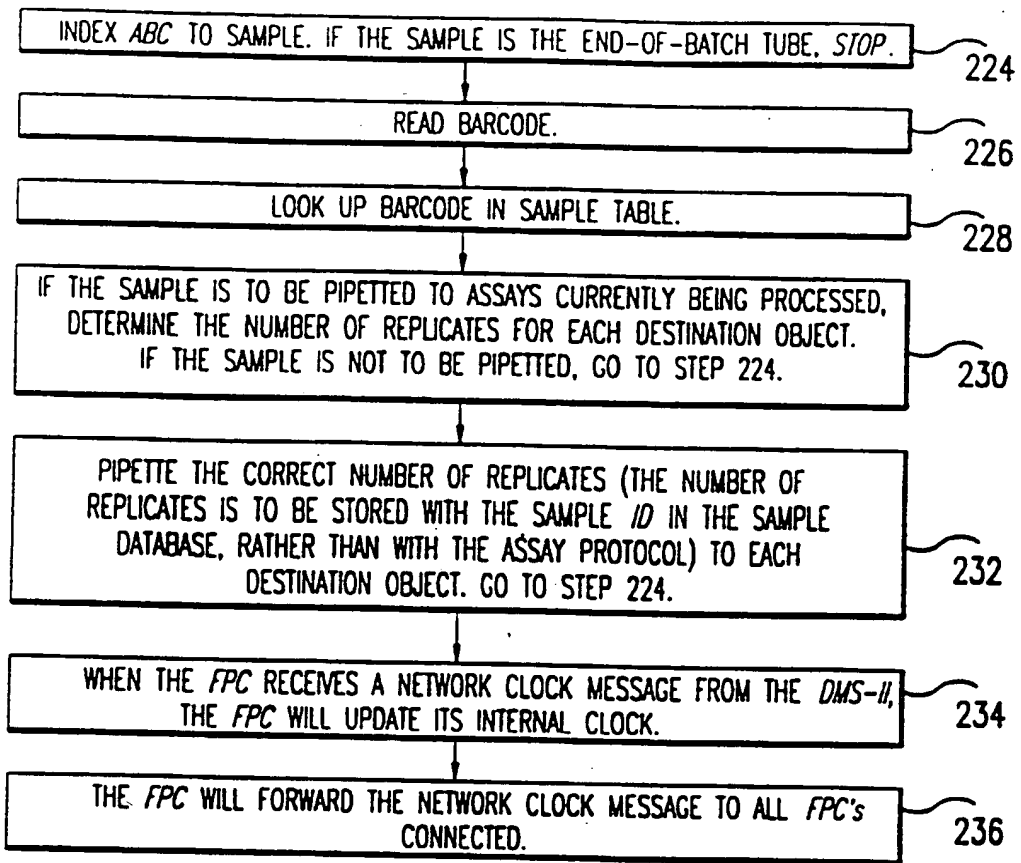


FIG.11

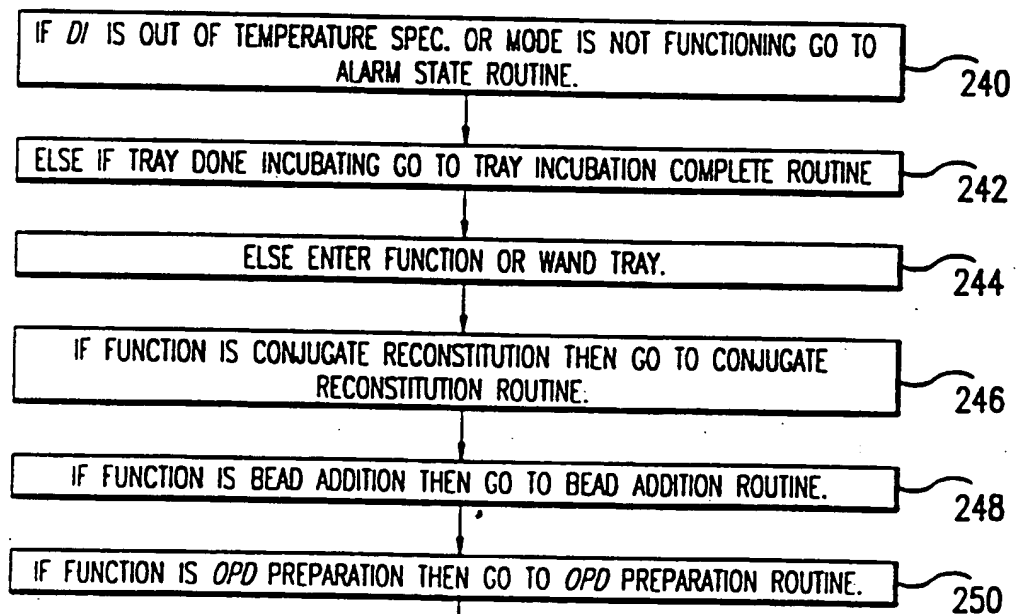


FIG.12

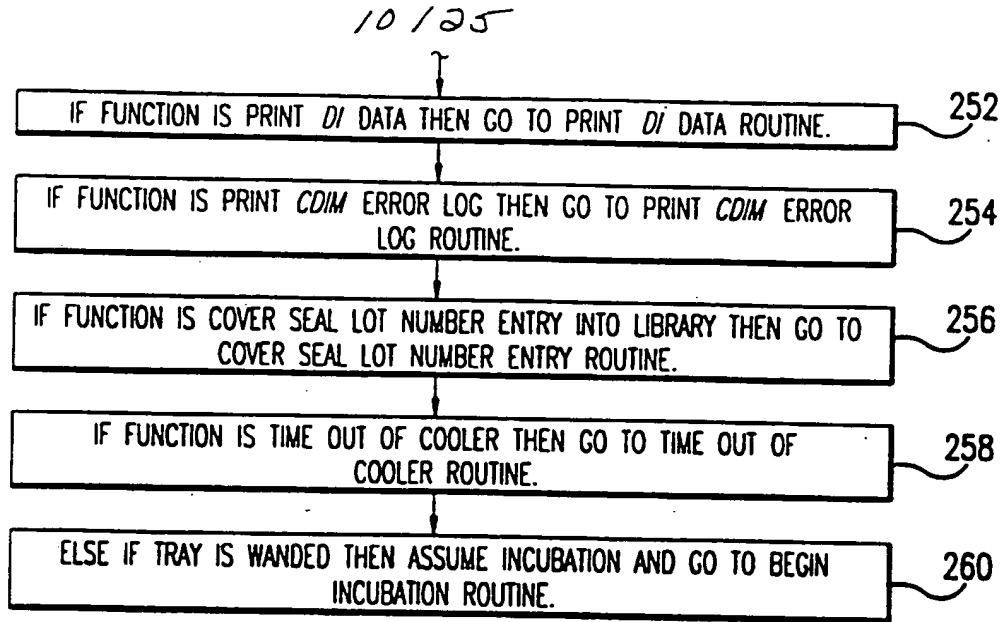


FIG.12 Cont.

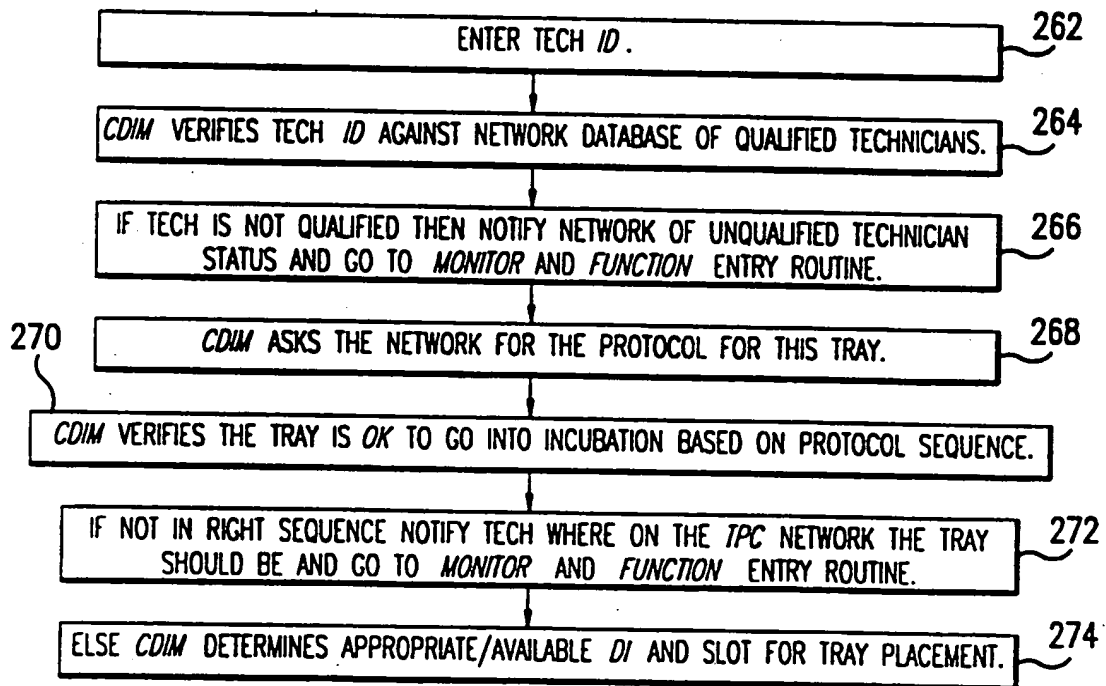


FIG.13a

11 / 25

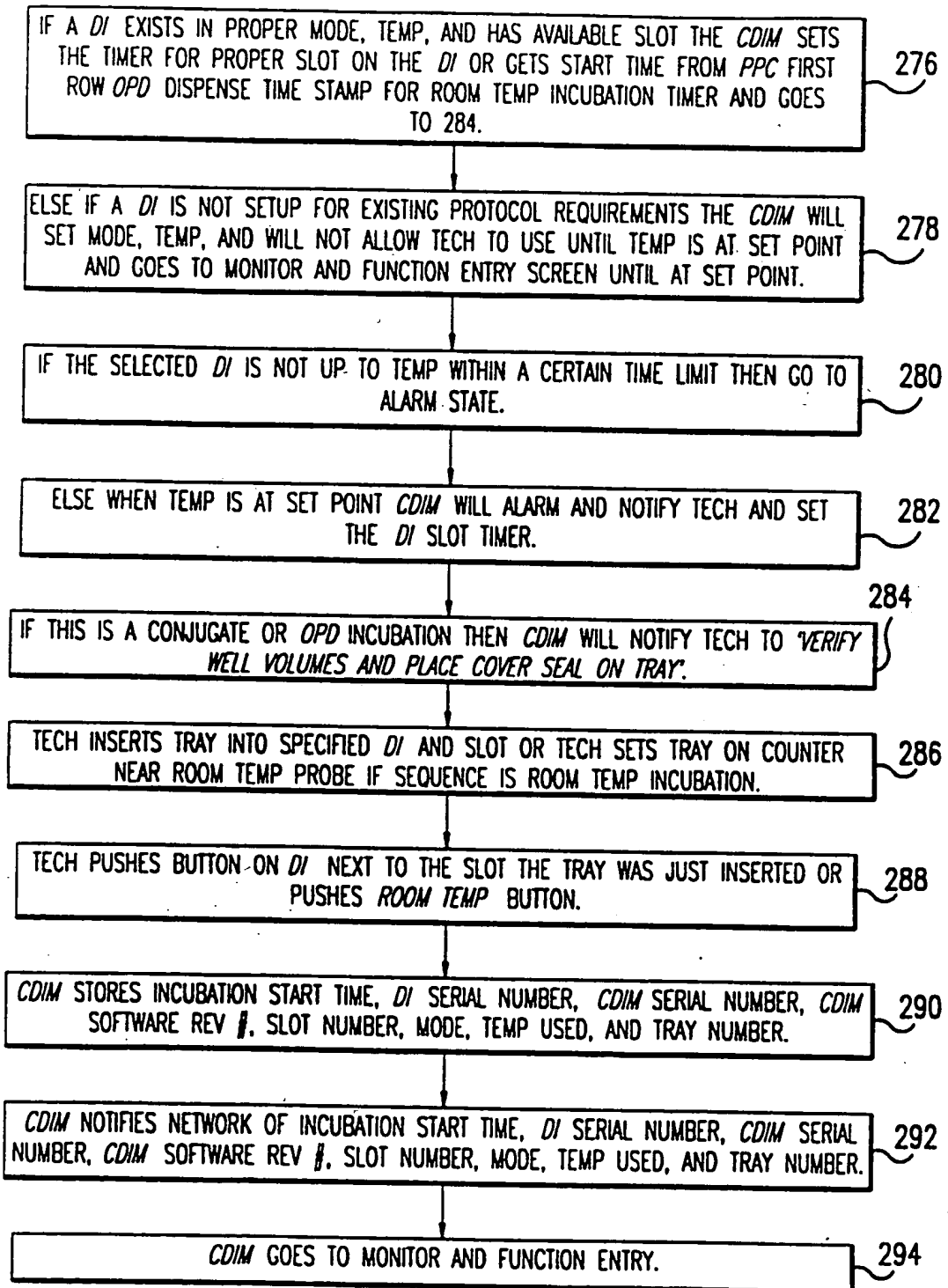


FIG.13b

12/25

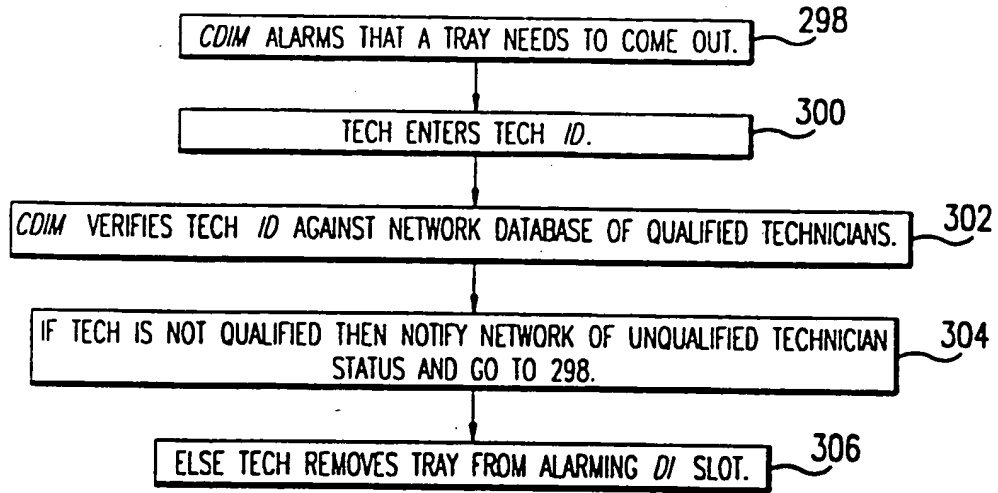


FIG. 14a

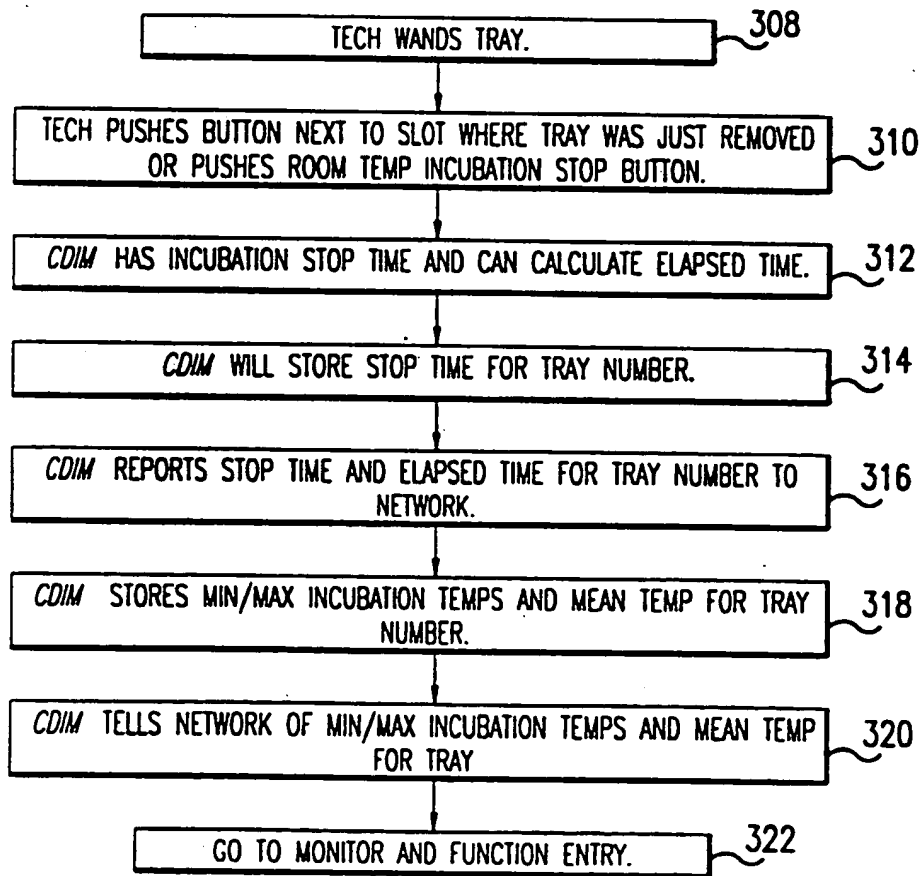


FIG. 14b

13/25

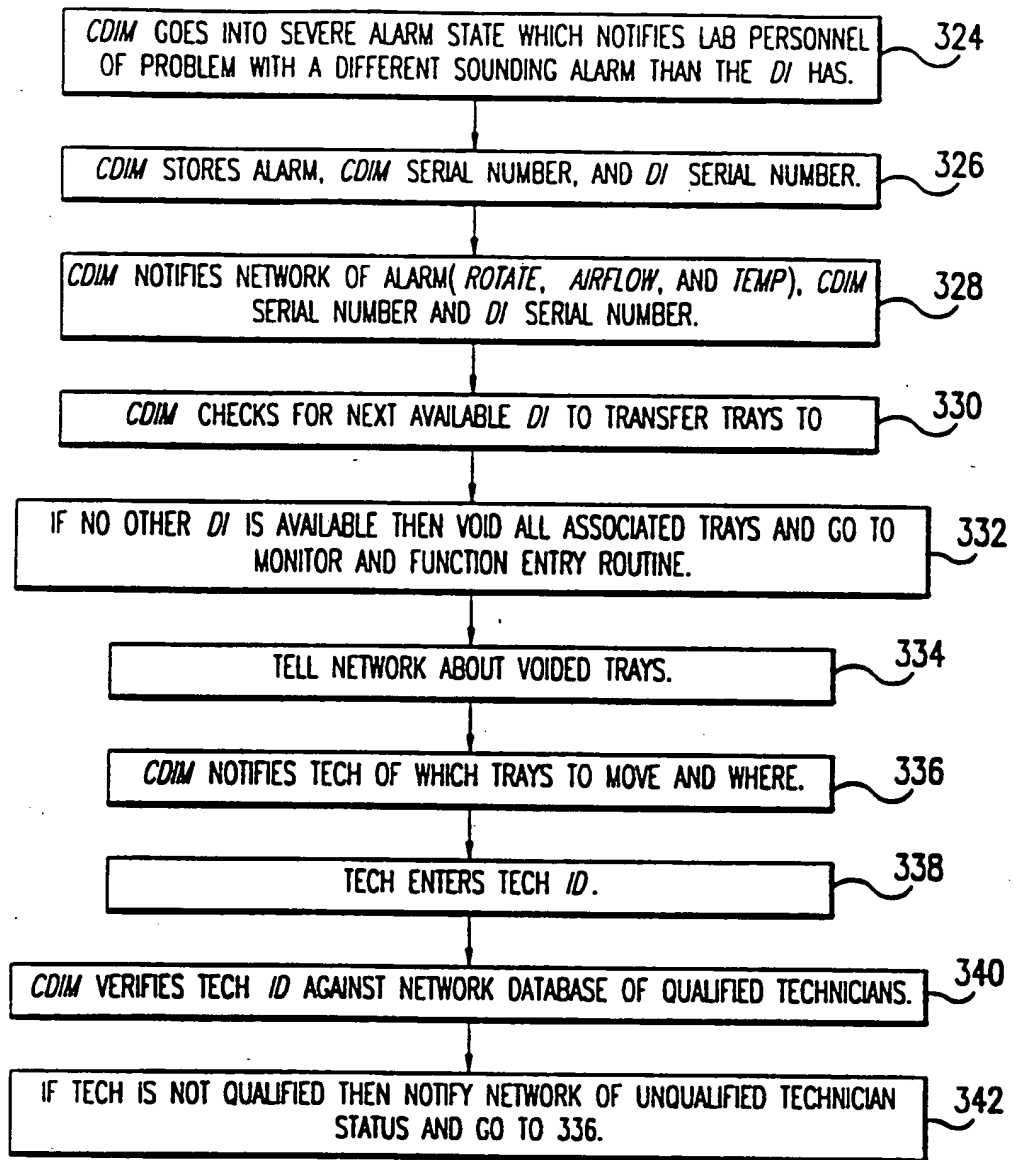


FIG. 15a

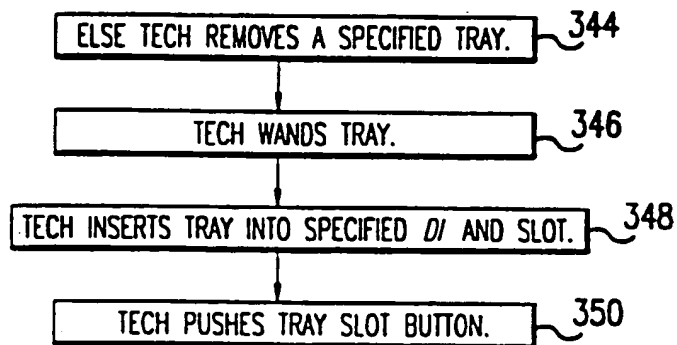


FIG. 15b

14/25

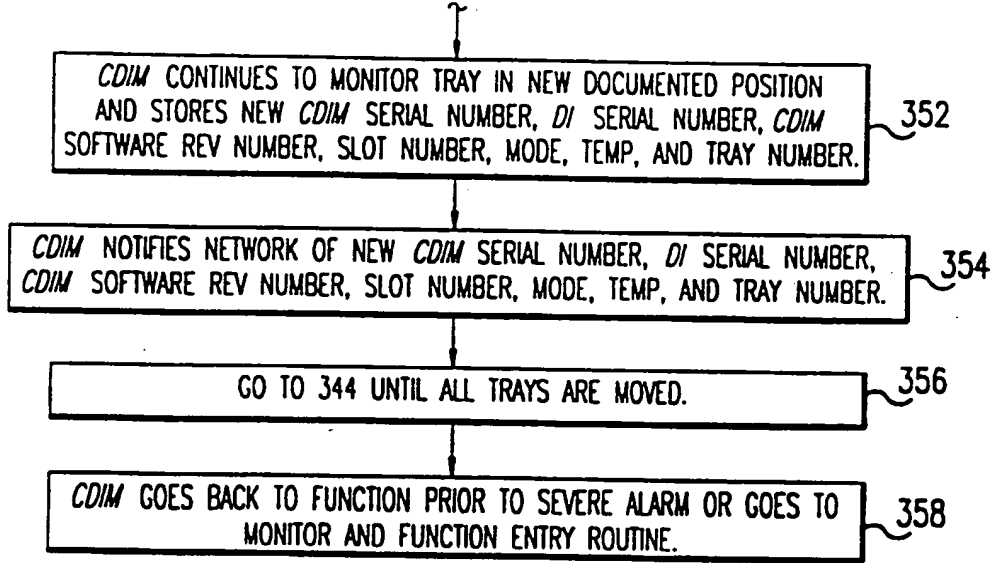


FIG.15b Cont.

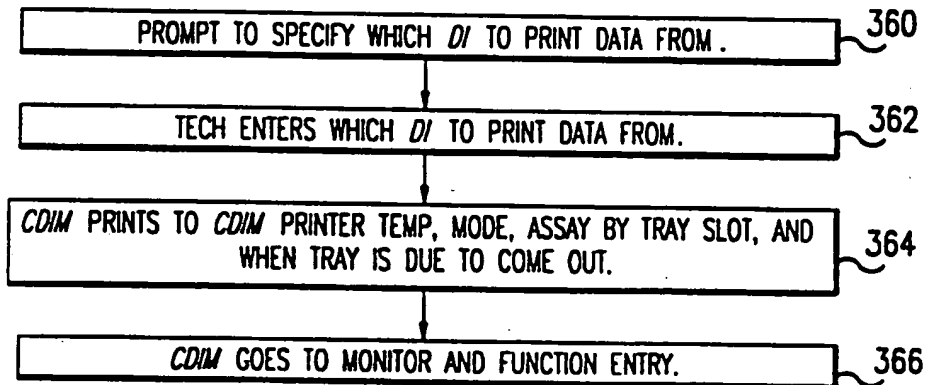


FIG.16

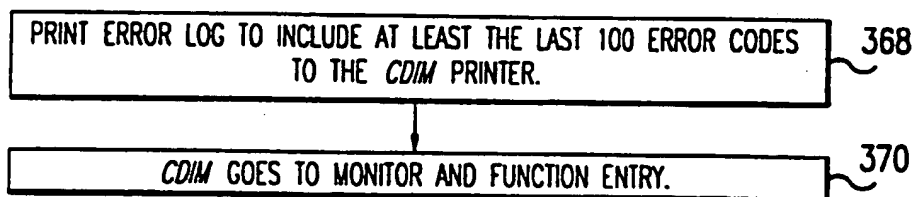


FIG.17



15/25

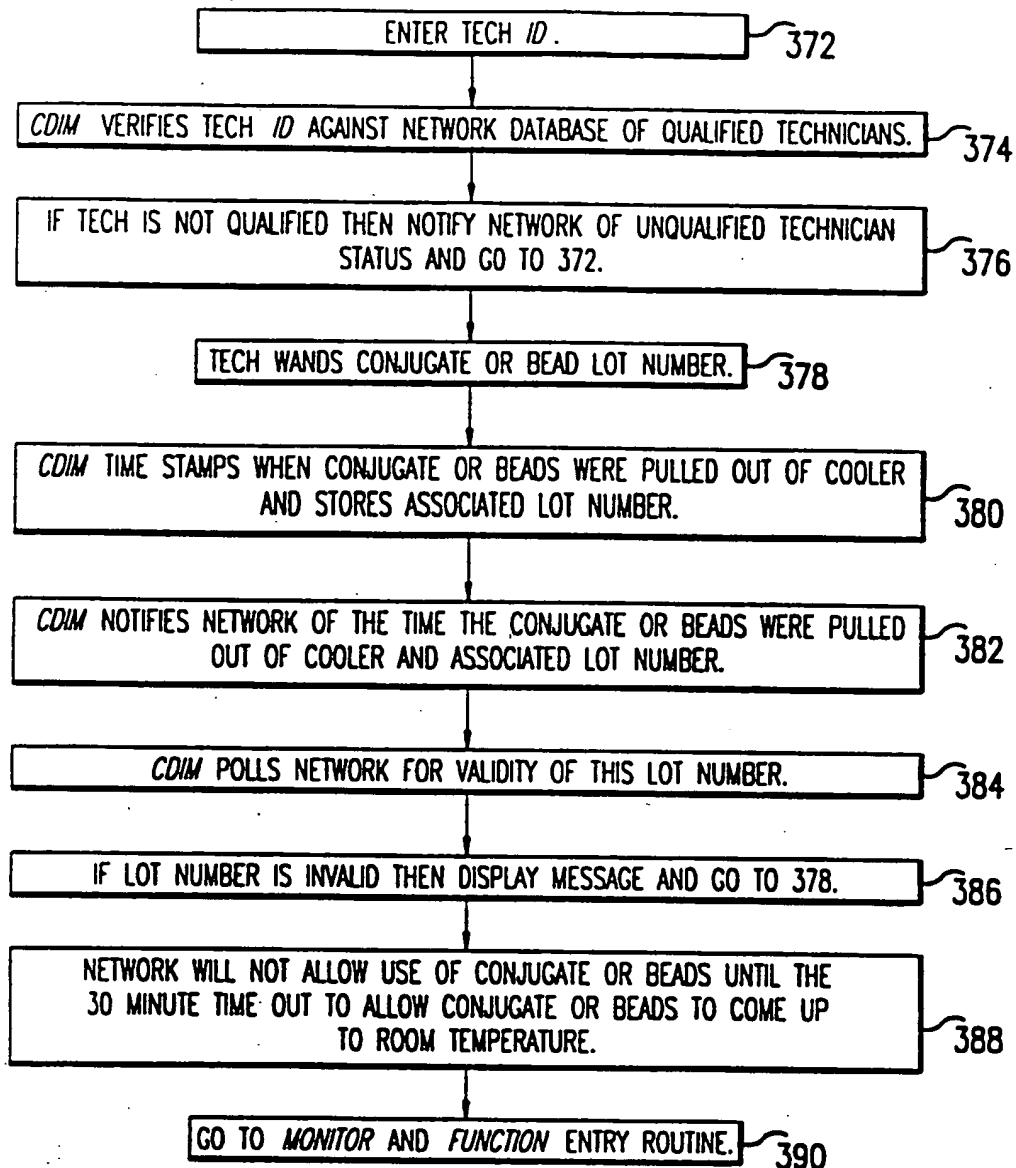


FIG.18

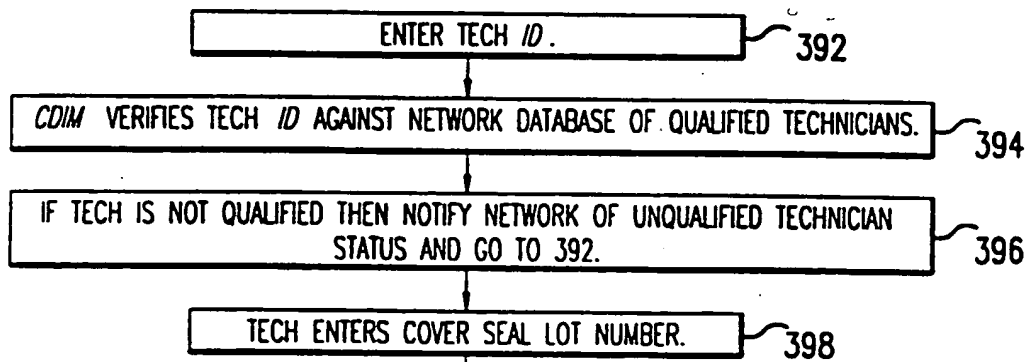


FIG.19

16/25

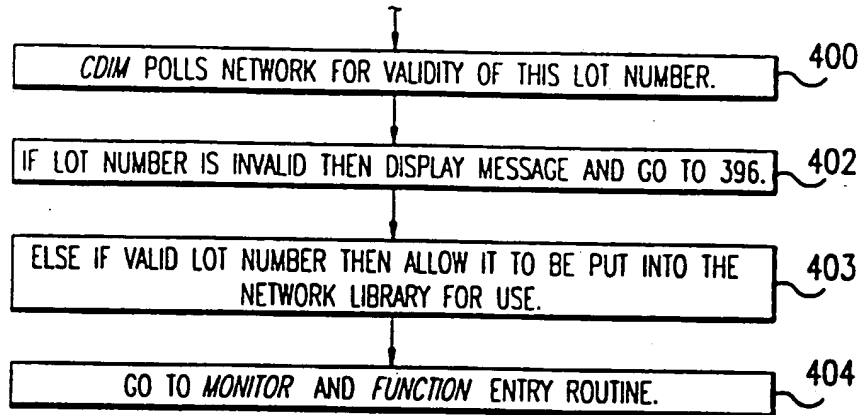


FIG.19 Cont.

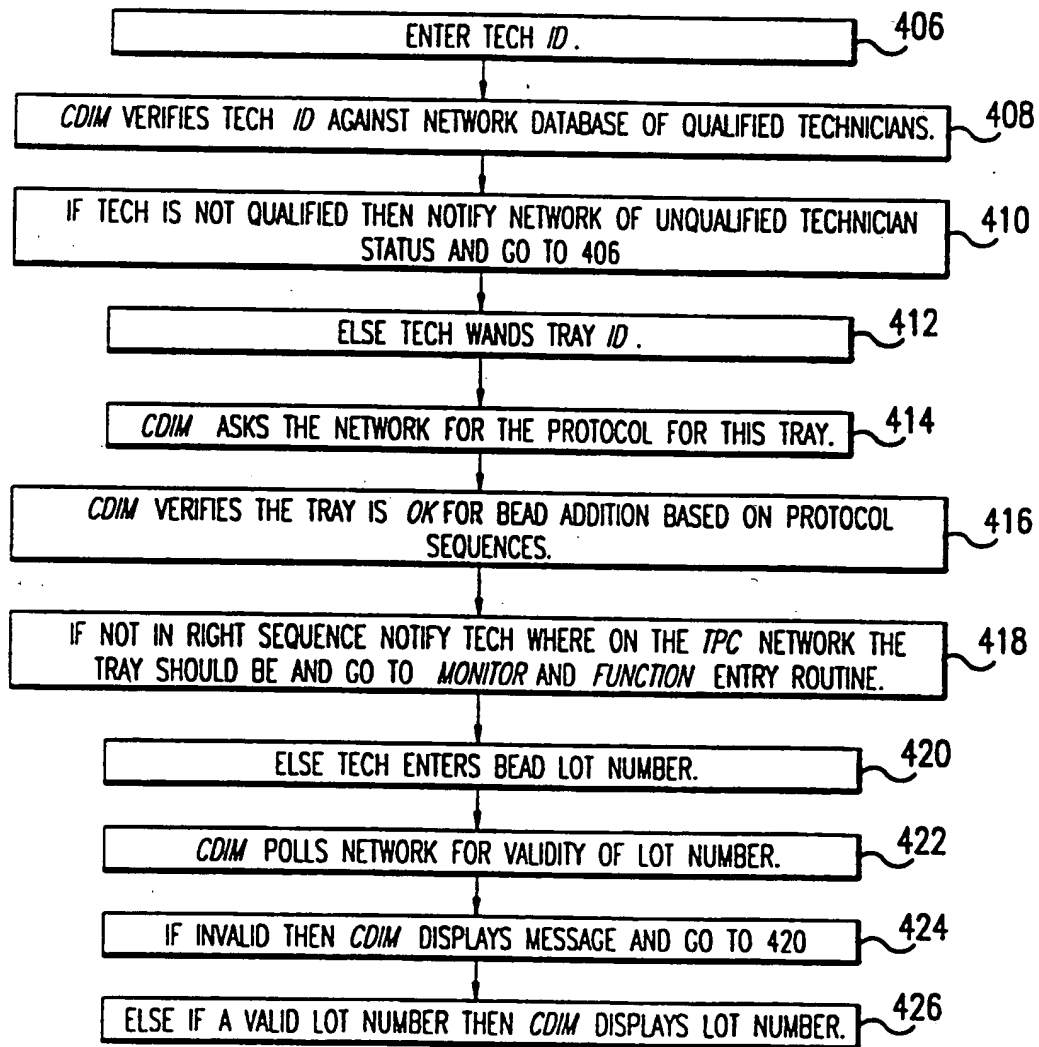


FIG.20a  
SUBSTITUTE SHEET (RULE 26)

17 / 25

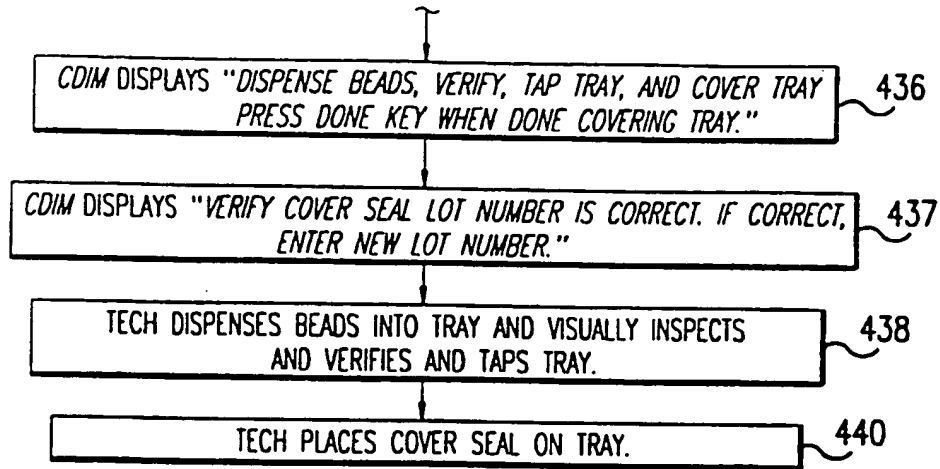


FIG.20a Cont.

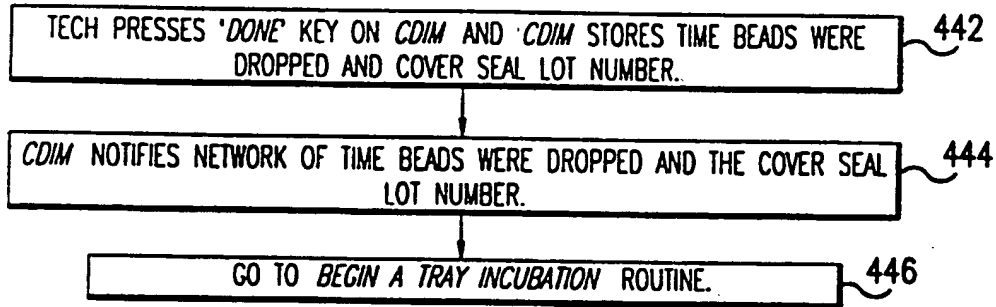


FIG.20b

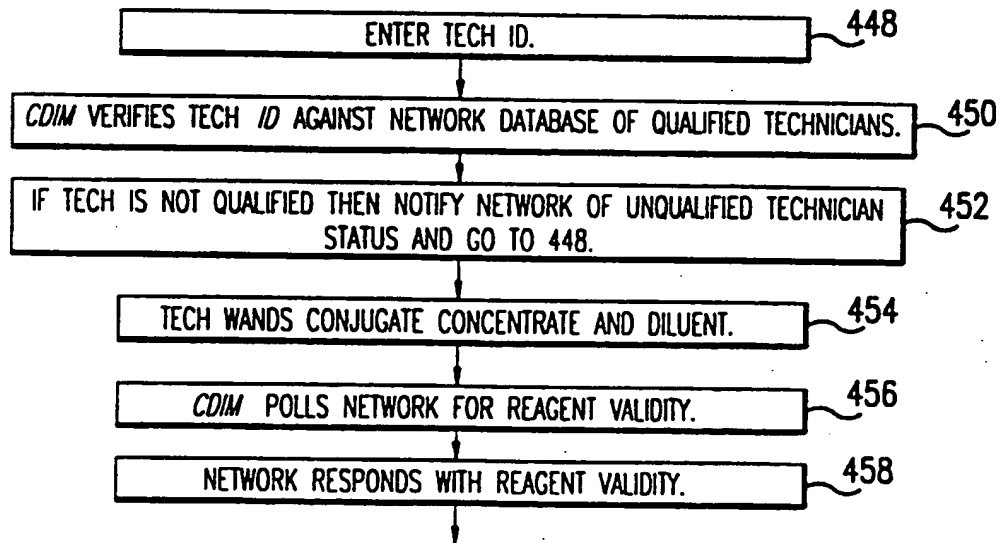


FIG.21a

18/25

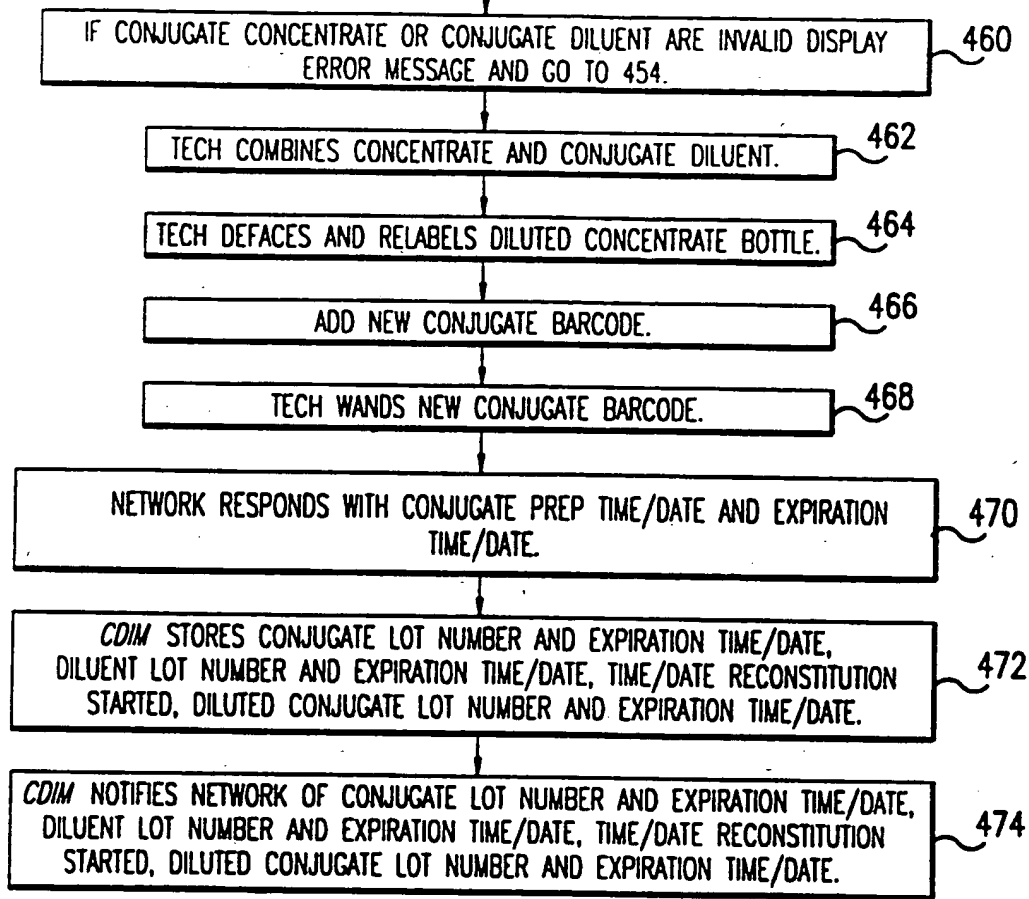


FIG.21a Cont.

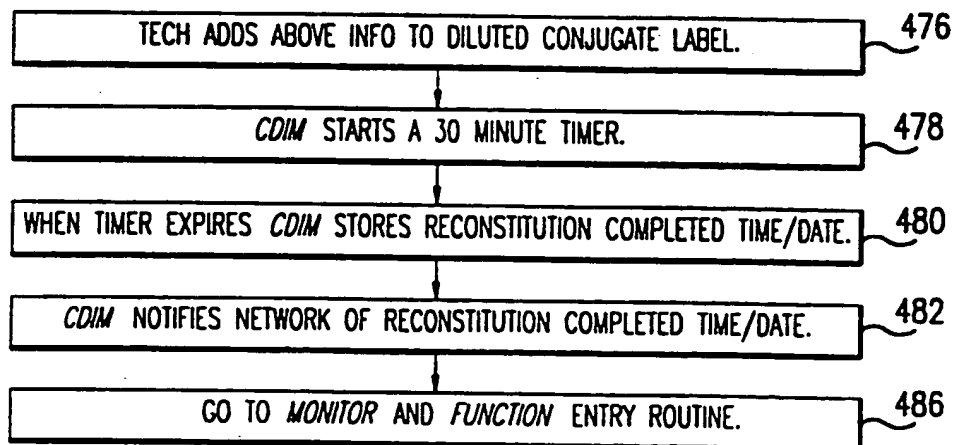


FIG.21b

19125

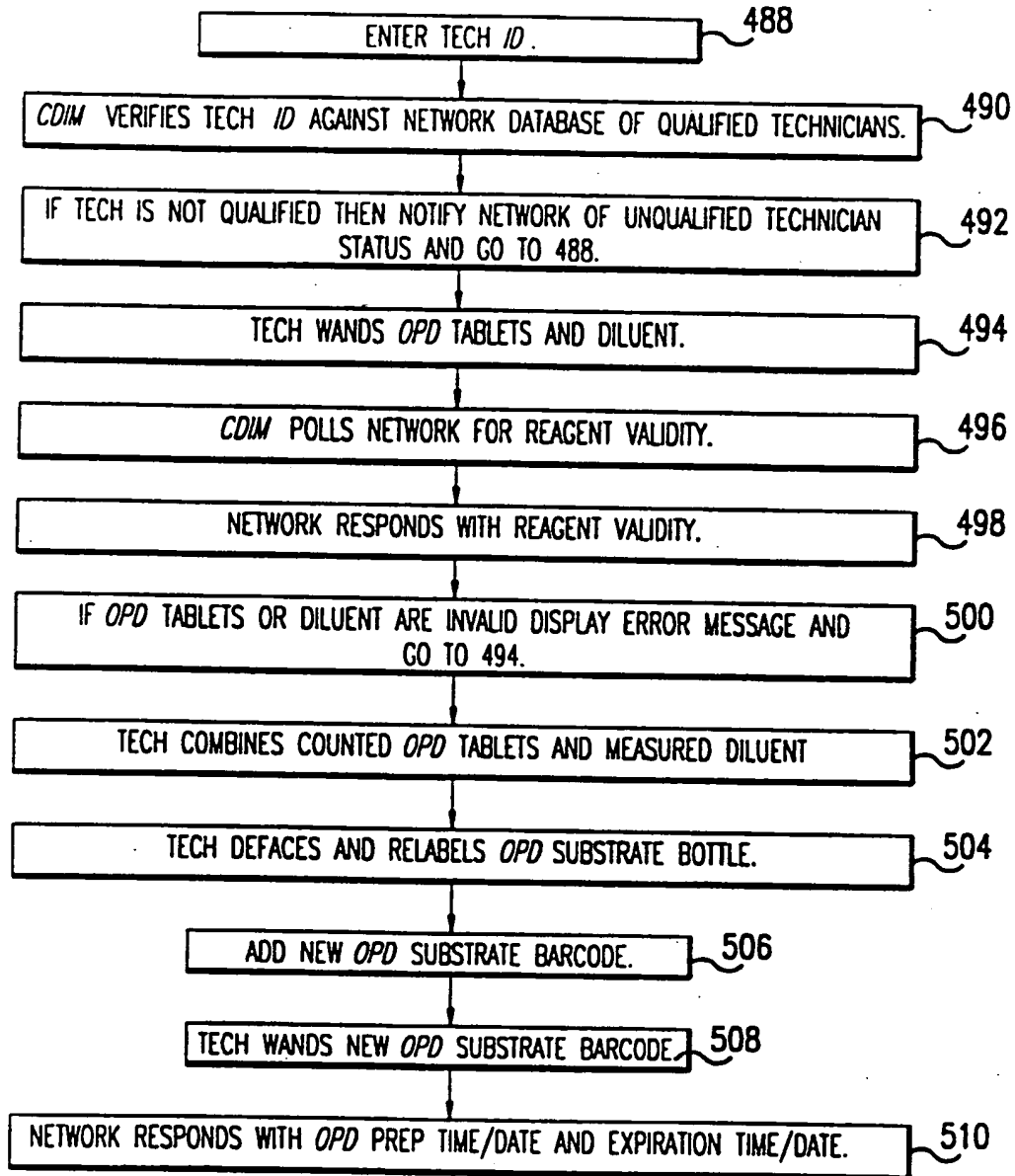


FIG.22a

20 / 25

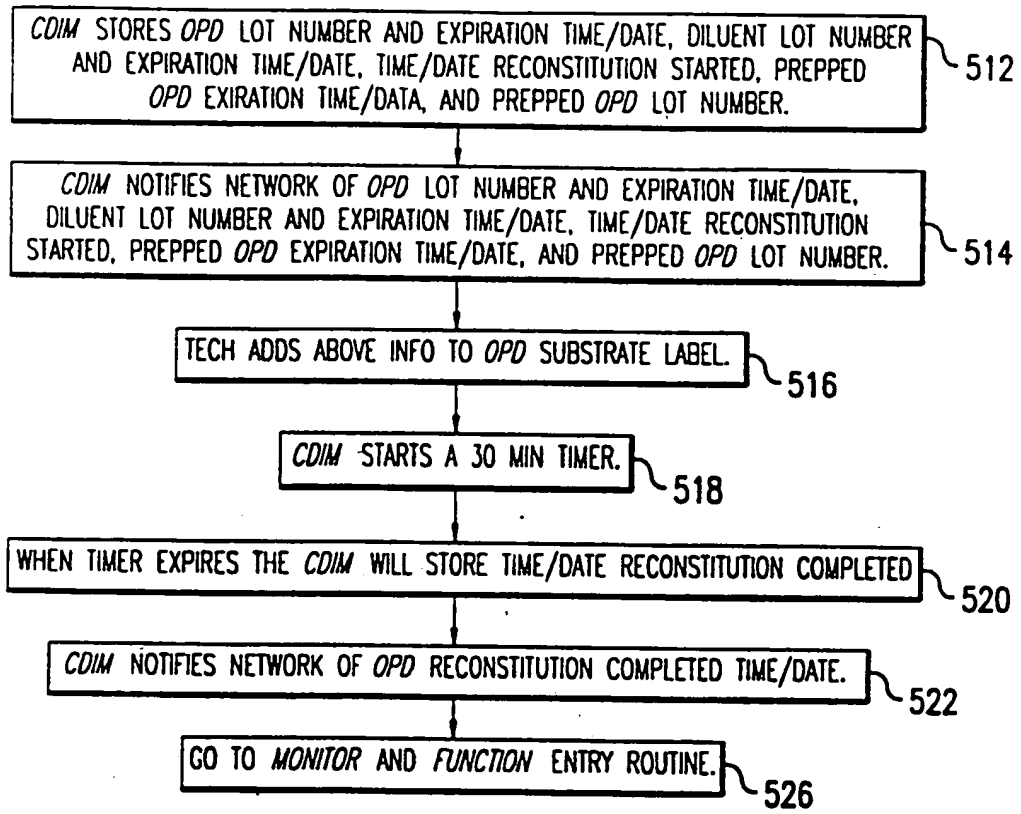


FIG. 22b

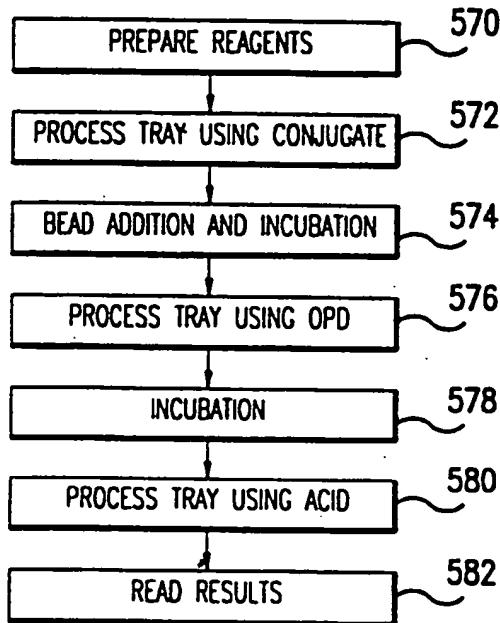


FIG. 23

21 / 25

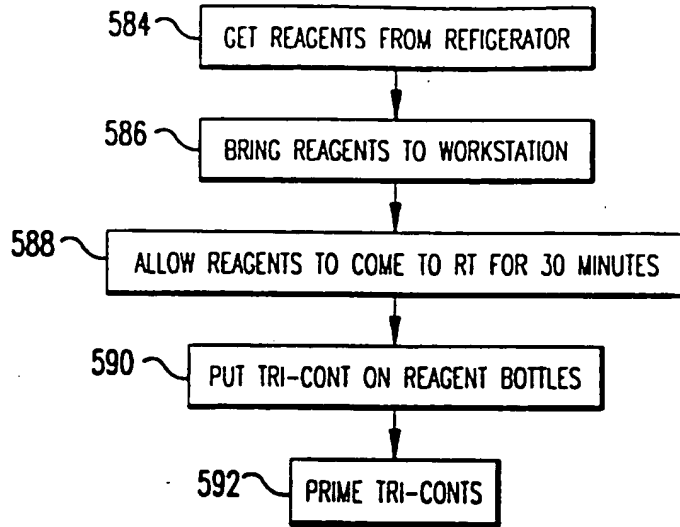


FIG.24

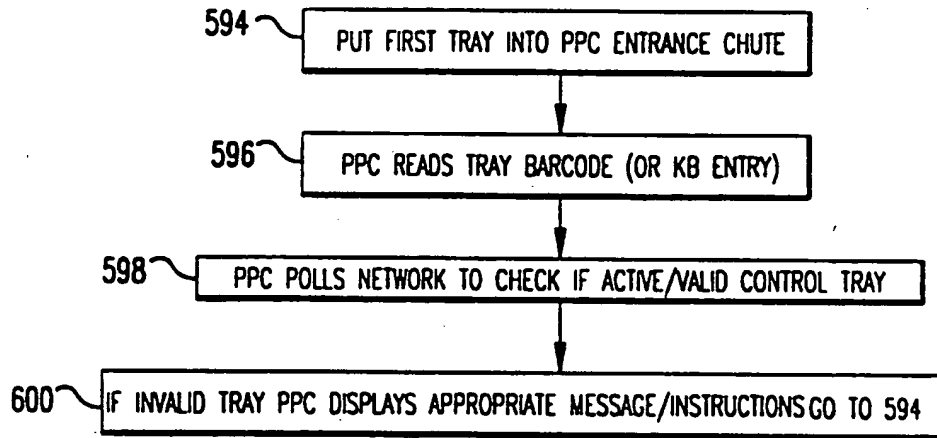


FIG.25

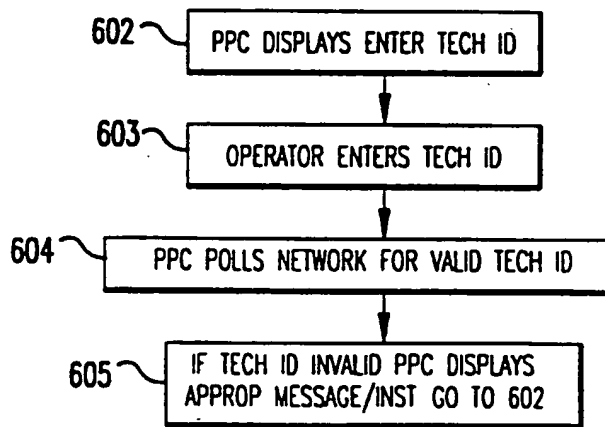


FIG.26

22 / 25

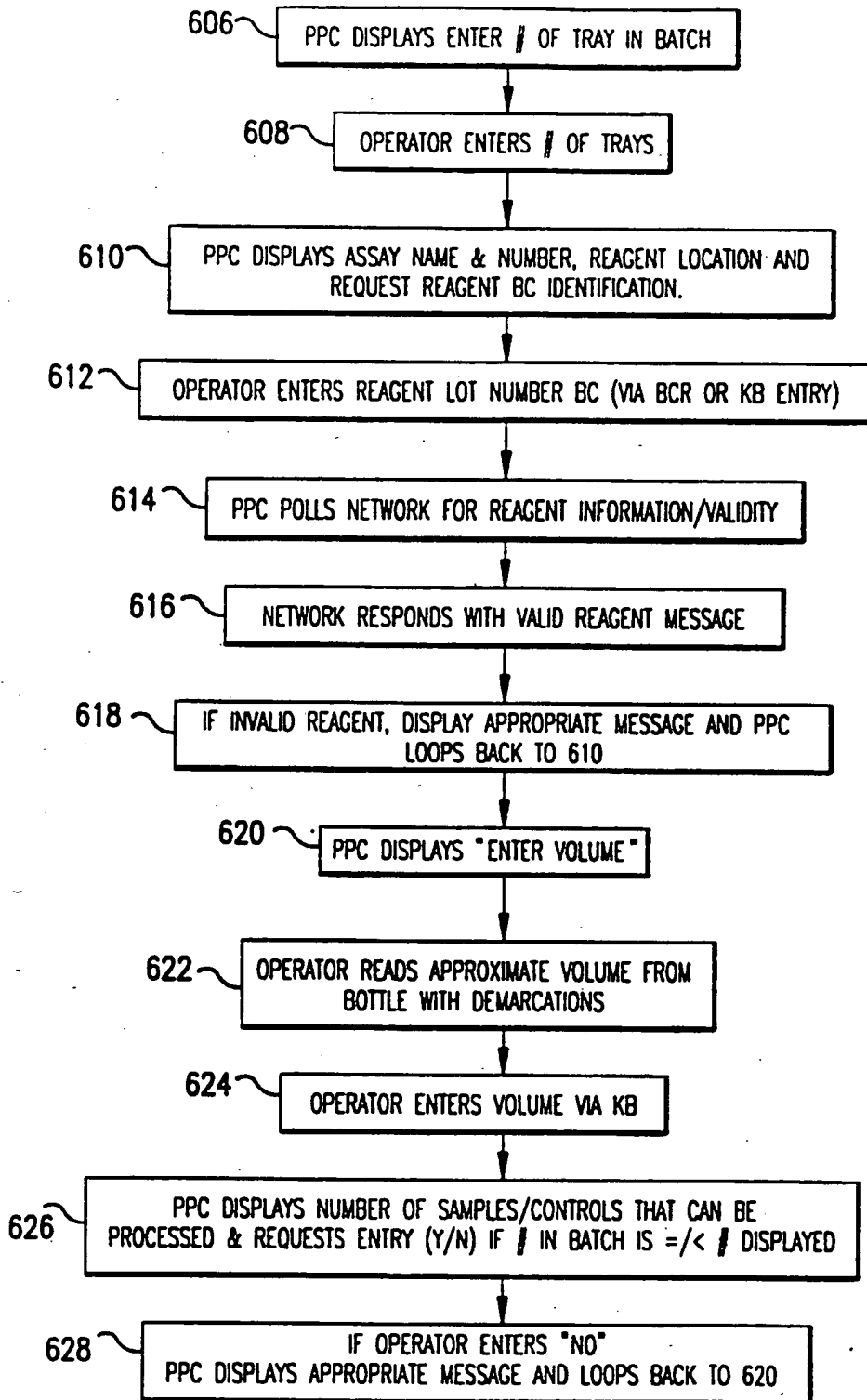


FIG.27a

SUBSTITUTE SHEET (RULE 26)



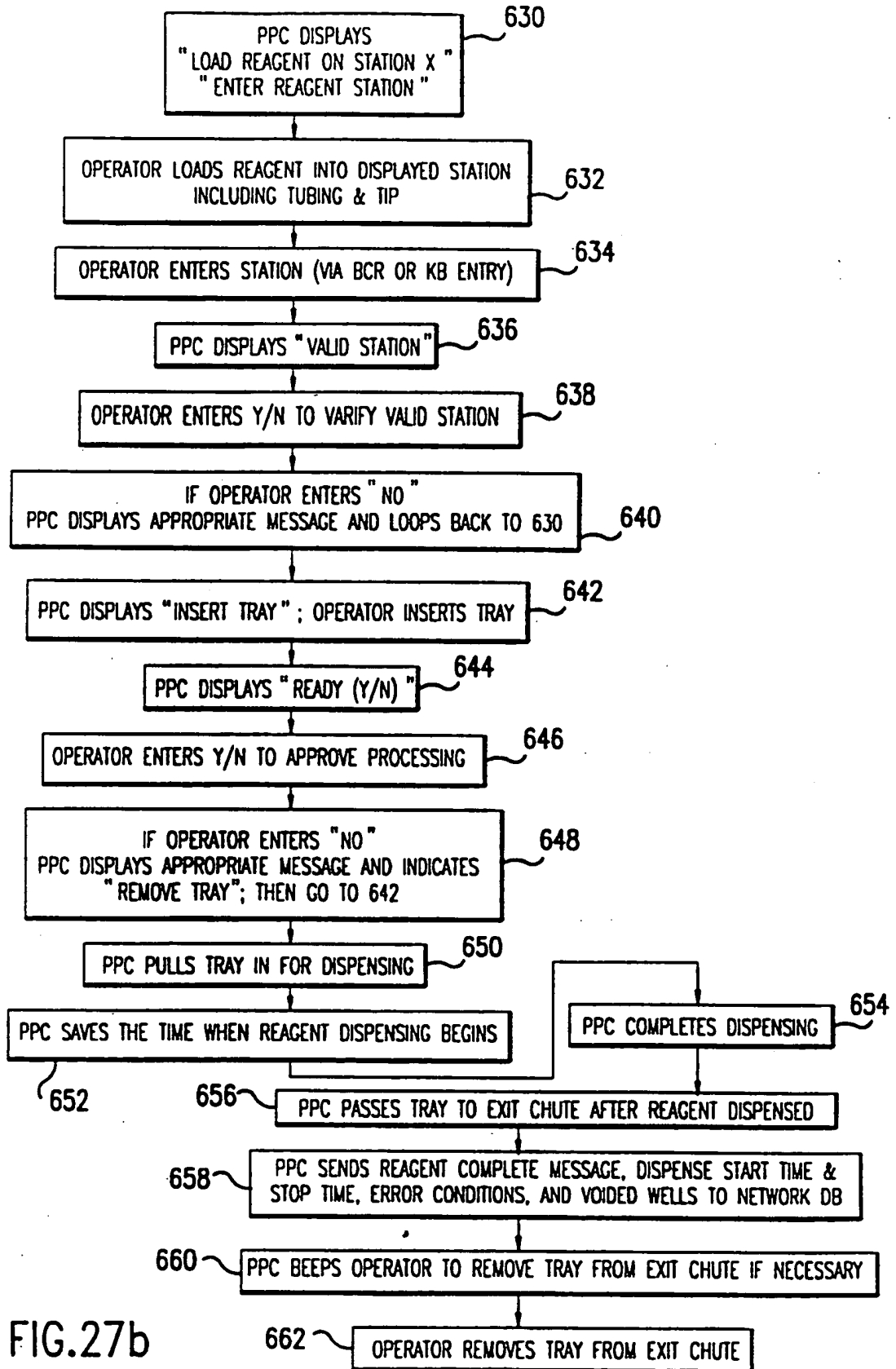


FIG.27b

24 / 25

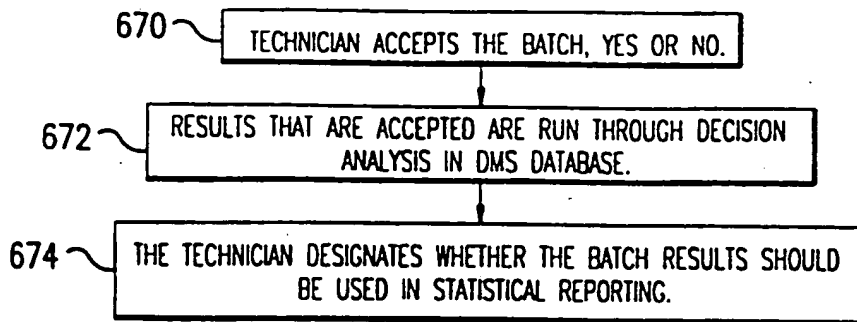


FIG.28

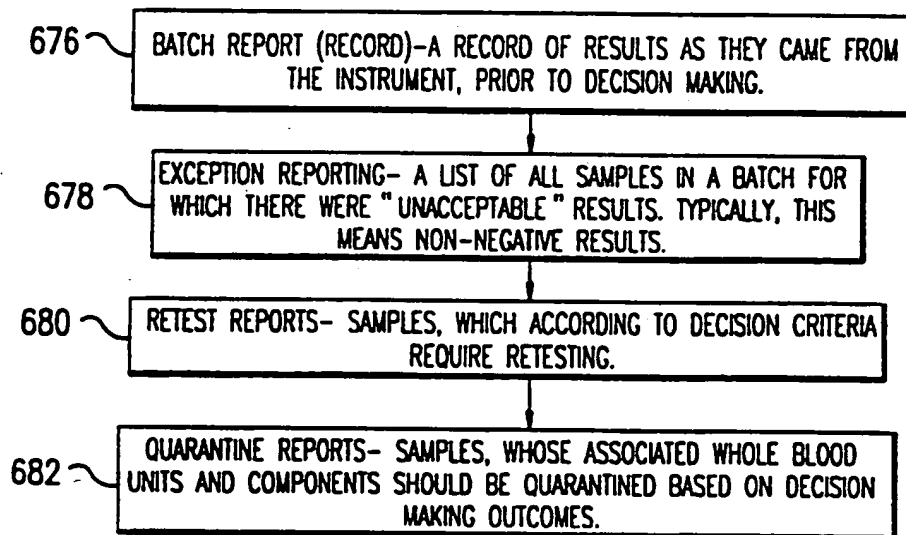


FIG.29

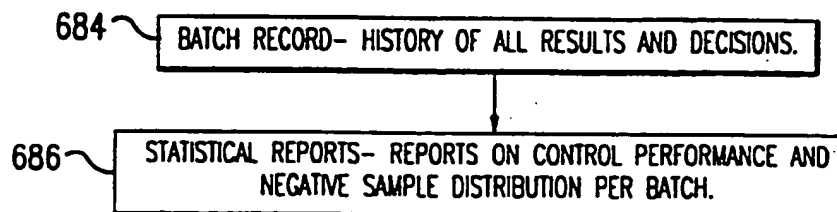


FIG.30

25/25

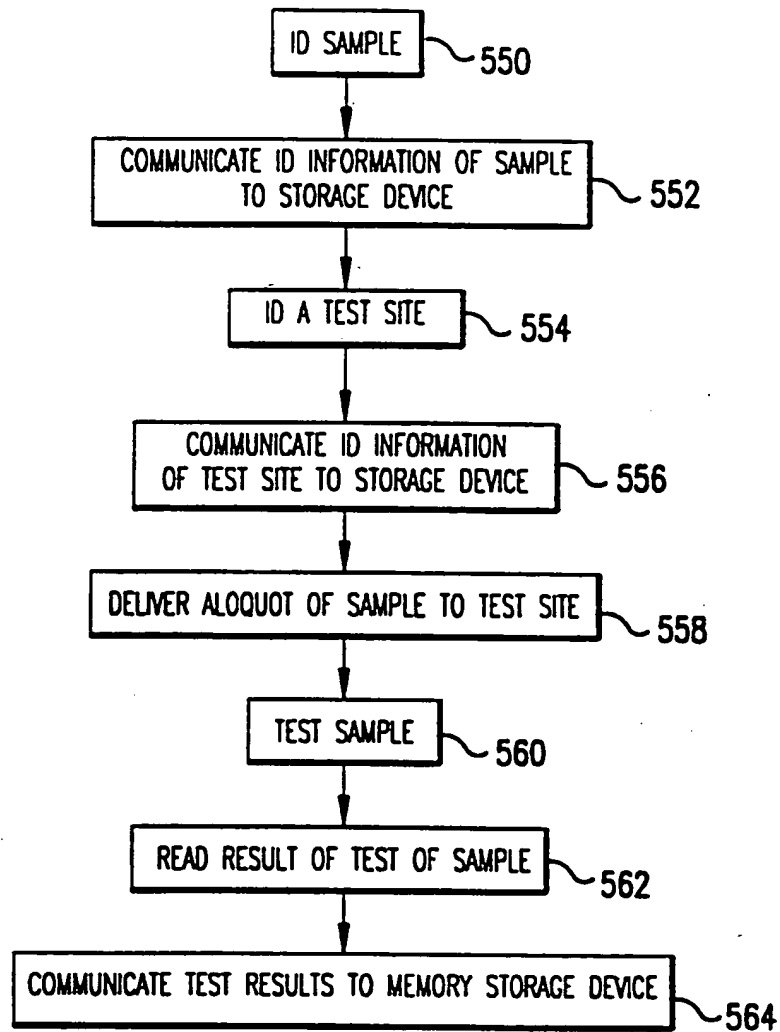


FIG.31

INTERNATIONAL SEARCH REPORT

International application No.  
PCT/US93/10525

<p><b>A. CLASSIFICATION OF SUBJECT MATTER</b>                  IPC(5) : G06F 45/42                  US CL : 364/413.01                  According to International Patent Classification (IPC) or to both national classification and IPC</p>																				
<p><b>B. FIELDS SEARCHED</b>                  Minimum documentation searched (classification system followed by classification symbols)                  U.S. : 364/413.01, 413.07, 413.08, 413.09, 413.10, 413.11, 468, 478, 497, 500                  Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched                  Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)</p>																				
<p><b>C. DOCUMENTS CONSIDERED TO BE RELEVANT</b></p> <table border="1"> <thead> <tr> <th>Category*</th> <th>Citation of document, with indication, where appropriate, of the relevant passages</th> <th>Relevant to claim No.</th> </tr> </thead> <tbody> <tr> <td>Y, P</td> <td>US, A, 5,231,585 (Kobayashi et al) 27 July 1993, abstract and figs. 1, 13-16B.</td> <td>1-7</td> </tr> <tr> <td>A, E</td> <td>US, A, 5,282,139 (Kobayashi) 25 January 1994, abstract and fig. 1.</td> <td>1-7</td> </tr> <tr> <td>A, E</td> <td>US, A, 5,282,149 (Grandone et al) 25 January 1994, abstract and fig. 1.</td> <td>1-7</td> </tr> <tr> <td>A</td> <td>US, A, 5,122,342 (McCulloch et al) 16 June 1992, abstract and figs. 1-2.</td> <td>1-7</td> </tr> <tr> <td>A, E</td> <td>US, A, 5,273,392 (Bernard, II et al) 28 December 1993, abstract and figs. 2 and 41-42.</td> <td>1-7</td> </tr> </tbody> </table>			Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.	Y, P	US, A, 5,231,585 (Kobayashi et al) 27 July 1993, abstract and figs. 1, 13-16B.	1-7	A, E	US, A, 5,282,139 (Kobayashi) 25 January 1994, abstract and fig. 1.	1-7	A, E	US, A, 5,282,149 (Grandone et al) 25 January 1994, abstract and fig. 1.	1-7	A	US, A, 5,122,342 (McCulloch et al) 16 June 1992, abstract and figs. 1-2.	1-7	A, E	US, A, 5,273,392 (Bernard, II et al) 28 December 1993, abstract and figs. 2 and 41-42.	1-7
Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.																		
Y, P	US, A, 5,231,585 (Kobayashi et al) 27 July 1993, abstract and figs. 1, 13-16B.	1-7																		
A, E	US, A, 5,282,139 (Kobayashi) 25 January 1994, abstract and fig. 1.	1-7																		
A, E	US, A, 5,282,149 (Grandone et al) 25 January 1994, abstract and fig. 1.	1-7																		
A	US, A, 5,122,342 (McCulloch et al) 16 June 1992, abstract and figs. 1-2.	1-7																		
A, E	US, A, 5,273,392 (Bernard, II et al) 28 December 1993, abstract and figs. 2 and 41-42.	1-7																		
<p><input checked="" type="checkbox"/> Further documents are listed in the continuation of Box C. <input type="checkbox"/> See patent family annex.</p>																				
<table border="0"> <tr> <td>* Special categories of cited documents:</td> <td>"T"</td> <td>later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention</td> </tr> <tr> <td>"A" document defining the general state of the art which is not considered to be part of particular relevance</td> <td>"X"</td> <td>document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone</td> </tr> <tr> <td>"E" earlier document published on or after the international filing date</td> <td>"Y"</td> <td>document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art</td> </tr> <tr> <td>"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)</td> <td>"&amp;"</td> <td>document member of the same patent family</td> </tr> <tr> <td>"O" document referring to an oral disclosure, use, exhibition or other means</td> <td></td> <td></td> </tr> <tr> <td>"P" document published prior to the international filing date but later than the priority date claimed</td> <td></td> <td></td> </tr> </table>			* Special categories of cited documents:	"T"	later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention	"A" document defining the general state of the art which is not considered to be part of particular relevance	"X"	document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone	"E" earlier document published on or after the international filing date	"Y"	document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art	"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)	"&"	document member of the same patent family	"O" document referring to an oral disclosure, use, exhibition or other means			"P" document published prior to the international filing date but later than the priority date claimed		
* Special categories of cited documents:	"T"	later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention																		
"A" document defining the general state of the art which is not considered to be part of particular relevance	"X"	document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone																		
"E" earlier document published on or after the international filing date	"Y"	document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art																		
"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)	"&"	document member of the same patent family																		
"O" document referring to an oral disclosure, use, exhibition or other means																				
"P" document published prior to the international filing date but later than the priority date claimed																				
<p>Date of the actual completion of the international search 13 APRIL 1994</p>		<p>Date of mailing of the international search report <b>25 APR 1994</b></p>																		
<p>Name and mailing address of the ISA/US                  Commissioner of Patents and Trademarks                  Box PCT                  Washington, D.C. 20231                  Facsimile No. N.A.</p>		<p>Authorized officer  <b>ROBERT WEINHARDT</b>  <i>B. Hardt</i>                  Telephone No. (703) 305-9780</p>																		

## INTERNATIONAL SEARCH REPORT

Inter. onal application No.  
PCT/US93/10525

C (Continuation). DOCUMENTS CONSIDERED TO BE RELEVANT		
Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
A	US, A, 4,857,716 (Gombrich et al) 15 August 1989, abstract, figs. 1, 14-15 and col. 8 line 8 to col. 9 line 40.	1-7
A,E	US, A, 5,260,872 (Copeland et al) 09 November 1993, abstract.	1-7
A, P	US, A, 5,197,011 (Biemans et al) 23 March 1993, abstract and fig. 1.	1-7